



**Scherzo SM-C18**  
**Scherzo SW-C18**

150 x 4.6 mm



Application

**Ibuprofen, Chlorpheniramine**

イブプロフェン, クロルフェニラミン

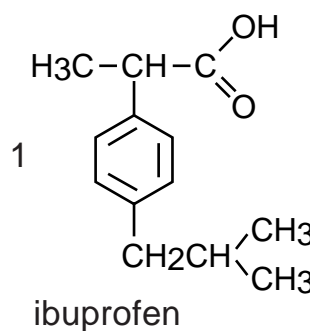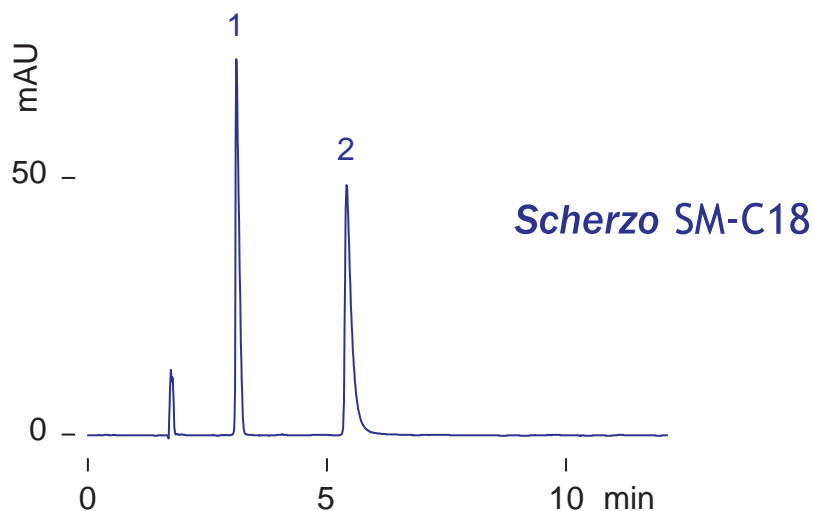

Scherzo SM-C18, 150 x 4.6 mm  
100mM ammonium acetate / methanol = 25 / 75  
1mL/min (12 MPa), 50 deg.C, 260 nm, 4 uL (1.3 - 8.7 ug)

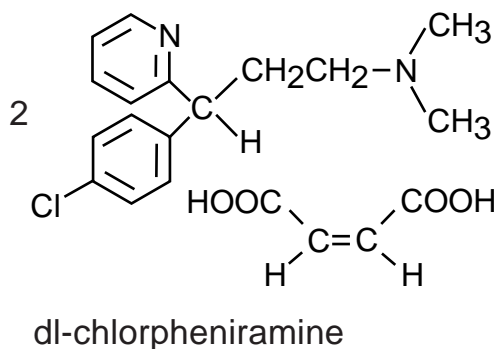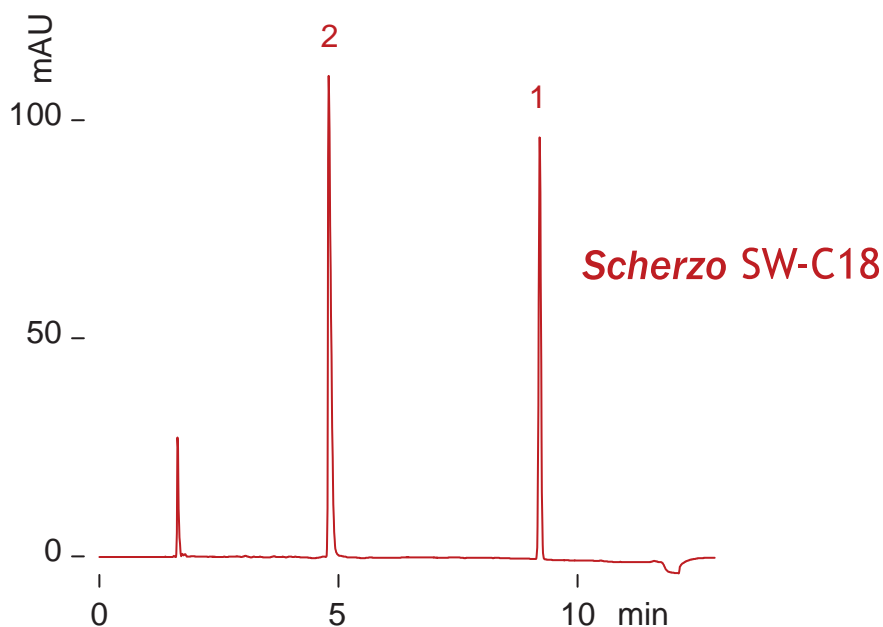

Scherzo SW-C18, 150 x 4.6 mm  
A: water / HCOOH = 100 / 0.3  
B: 200mM ammonium formate / acetonitrile = 25 / 75  
40-100 %B (0-10 min)  
1mL/min (9 MPa), 50 deg.C, 260 nm, 4 uL ( 1.3 - 8.7 ug)