

Cadenza CD-C18


100 x 2 mm

Application

LC-MS Application for Capsaicins

LC-MSアプリケーション（唐辛子成分、カプサイシン）

API4000 (TAKARA BIO INC.)

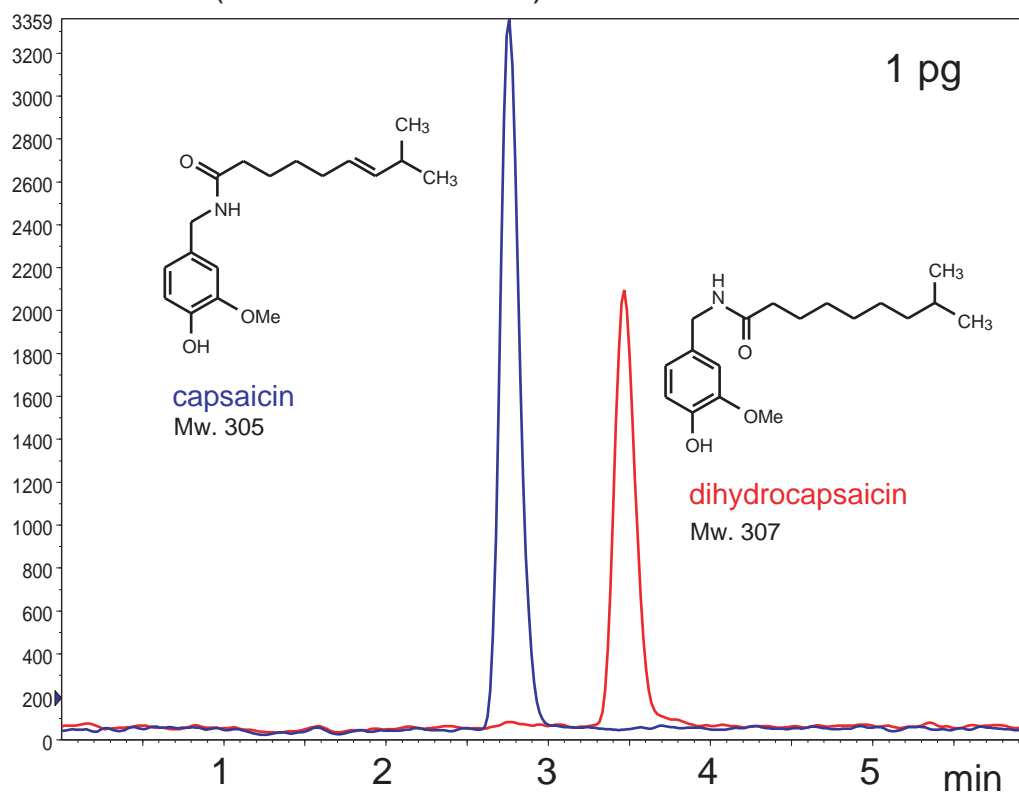

Cadenza CD-C18, 100 x 2 mm

0.1% acetic acid / 75% methanol, 0.2 mL/min, 40 °C, 5 μ L (200fg/ μ L)

API4000: ESI, MRM Positive

Q1/Q3:

capsaicin 306/137

dihydrocapsaicin 308/137

Courtesy of J.Watanabe, TAKARA BIO INC.