


Unison UK-Phenyl  
 Unison UK-C8  
 Unison US-C18  
 Unison UK-C18

250 x 4.6 mm

Application

Organic Acids  
 有機酸

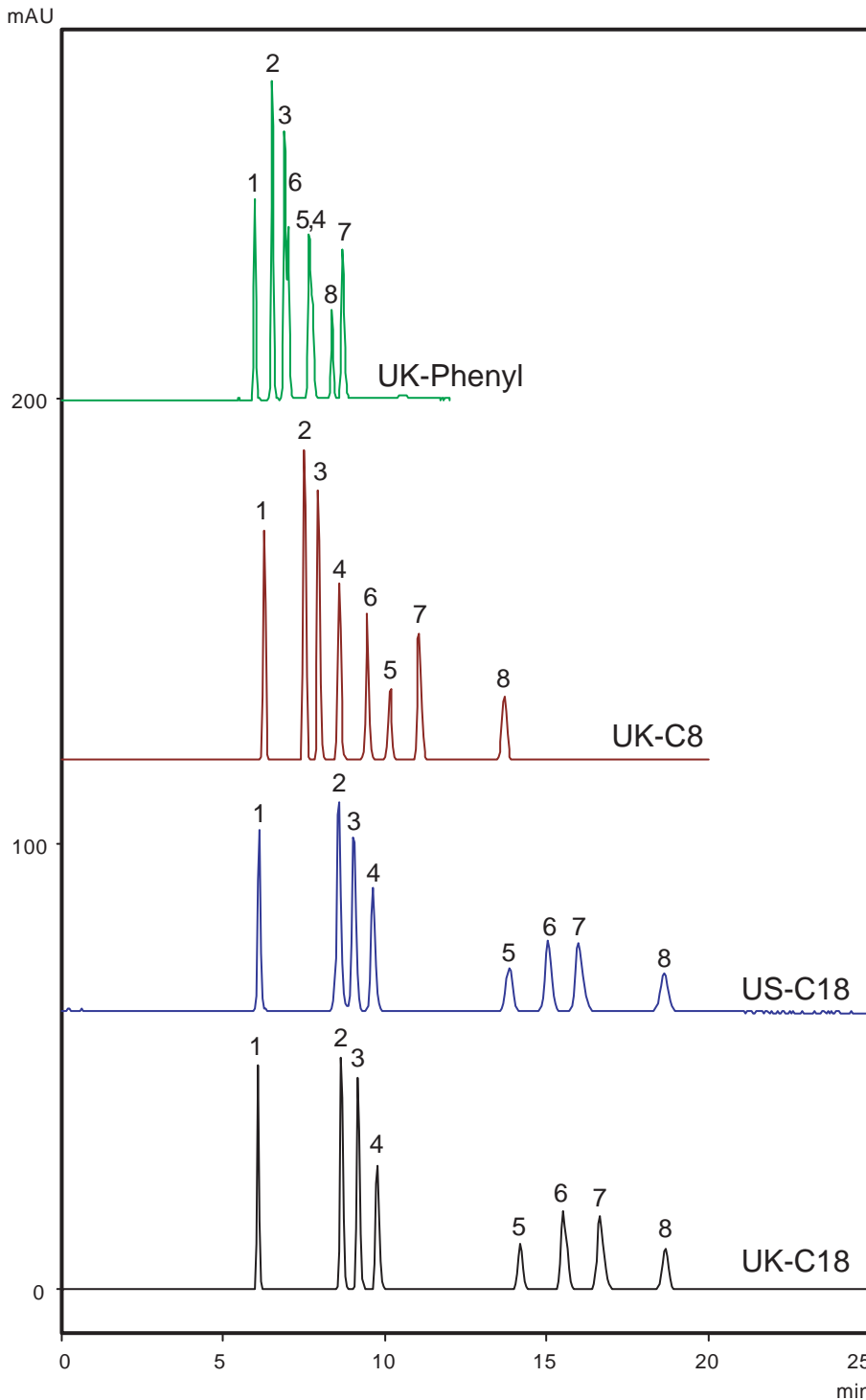

- 1 HCOOH  
formic acid
- 2 HOOCCH<sub>2</sub>COOH  
malonic acid
- 3 CC(O)C(=O)O  
lactic acid
- 4 CH<sub>3</sub>COOH  
acetic acid
- 5 OC(=O)C=CC(=O)O  
maleic acid
- 6 OC(=O)C(O)C(=O)O  
citric acid
- 7 HOOC(CH<sub>2</sub>)<sub>2</sub>COOH  
succinic acid
- 8 OC(=O)C=CC(=O)O  
fumaric acid

Unison, 250 x 4.6 mm  
 20 mM H<sub>3</sub>PO<sub>4</sub>  
 0.6 mL/min, 37 °C  
 UV at 210 nm