

Unison UK-C8



150 x 4.6 mm

Application

**Fatty Acid ADAM Derivatives**

脂肪酸のADAM誘導体

ADAM: 9-anthryldiazomethane

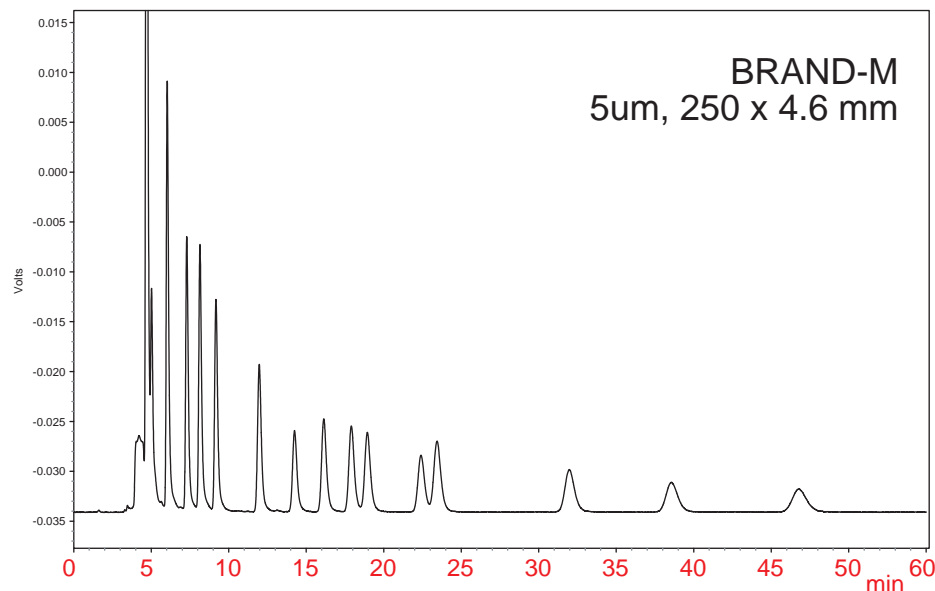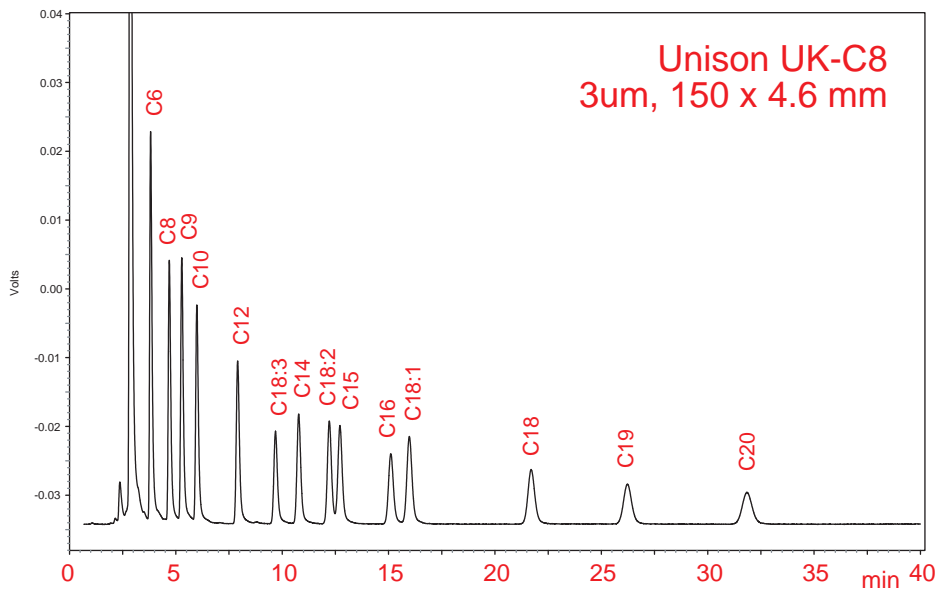

- C6** caproic acid  
n-hexanoic acid
- C8** caprylic acid  
n-octanoic acid
- C9** pelargonic acid  
n-nonanoic acid
- C10** capric acid  
n-decanoic acid
- C12** lauric acid  
n-dodecanoic acid
- C18:3** alpha-linolenic acid  
cis-9,12,15-octadecatrienoic acid
- C14** myristic acid  
n-tetradecanoic acid
- C18:2** linoleic acid  
cis-9,12-octadecadienoic acid
- C15** pentadecanoic acid  
n-pentadecanoic aci
- C16** palmitic acid  
n-hexadecanoic acid
- C18:1** oleic acid  
cis-9-octadecenoic acid
- C18** stearic acid  
n-octadecanoic acid
- C19** nonadecanoic acid  
n-nonadecanoic acid
- C20** arachidic acid  
eicosanoic acid

Unison UK-C8, 150 x 4.6 mm  
water / acetonitrile = 15 / 85, 1.0 mL/min, 30 °C  
FLS Ex365 nm, Em412 nm