


Unison UK-C18

250 x 4.6 mm

Application

GDP-Sugars

GDP糖

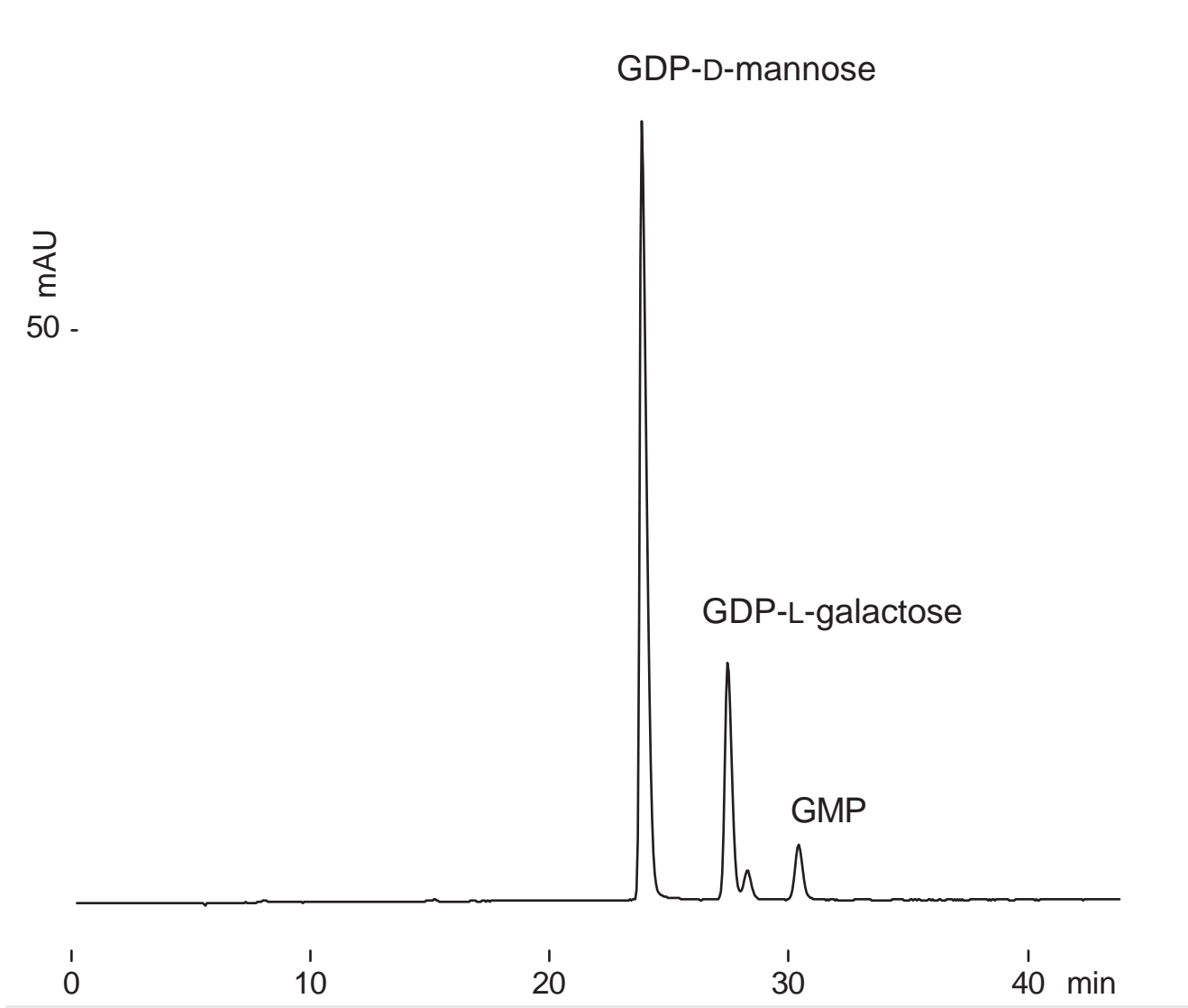

Unison UK-C18, 250 x 4.6 mm
 water / triethylamine / acetic acid = 100 / 0.2 / 0.1
 0.8 mL/min (13MPa), 37 deg.C, 260 nm, 5 uL

Ref) *Characterization of a GDP-D-mannose 3",5"-epimerase from rice*
 K. Watanabe, K. Suzuki, S. Kitamura, *Phytochemistry* 67 (2006) 338-346

Courtesy of Prof. Shinichi Kitamura, Osaka Pref. Univ., Japan