

Cadenza CD-C18



75 x 4.6 mm

Application


ABEE-Converted Oligosaccharide

ABEE標識化オリゴ糖

< α -Linked glucooligosaccharide >

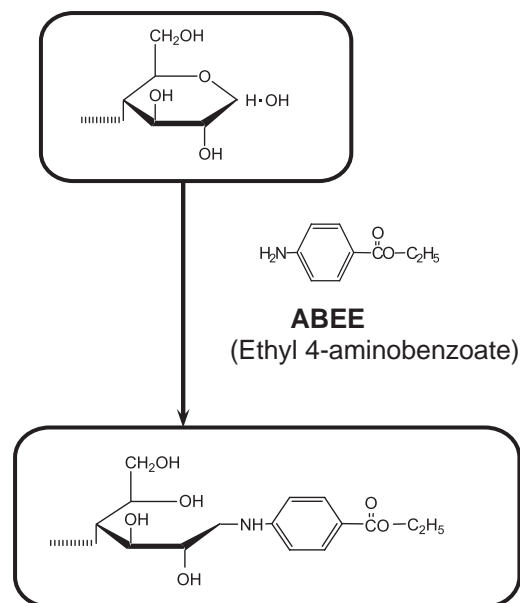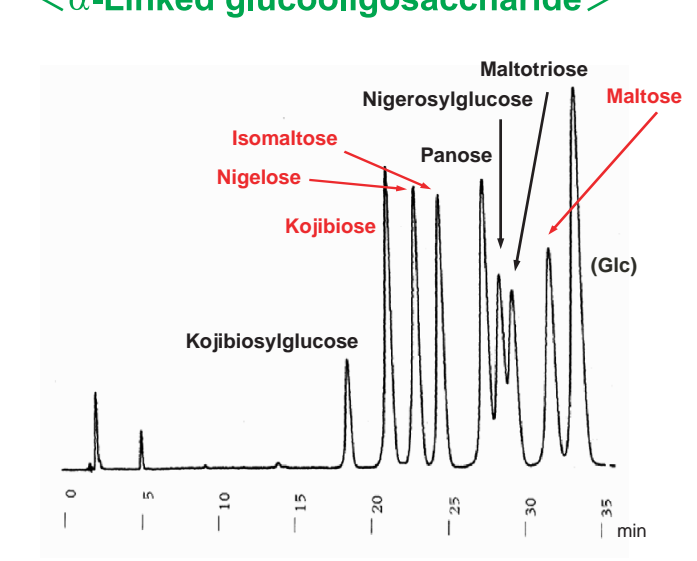

< α -Galactobiose (α -Gal2) >

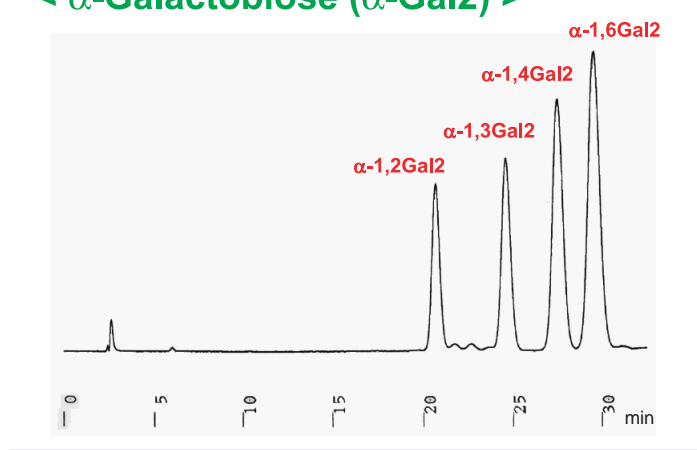

Cadenza CD-C18, 75 x 4.6 mm

0.1M ammonium acetate buffer (pH 4.0) / acetonitrile = 89 / 11 (v/v)

0.5 mL/min, 20 deg.C, 305 nm

Courtesy of Dr. Akiko Yamashita and Dr. Hiroyuki Hashimoto, Shinshu Univ., Japan