


Unison UK-Phenyl
 Unison UK-C8
 Cadenza CD-C18
 Cadenza CL-C18
 Unison UK-C18

75 x 3 mm


Application

Veterinary Medicine (Synthetic Antimicrobial)

動物用医薬品(合成抗菌剤)

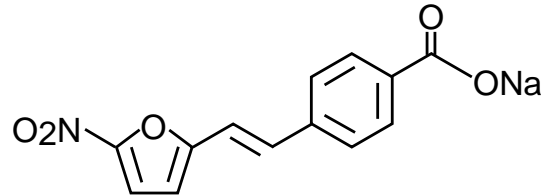

nifurstyreneate monosodium salt

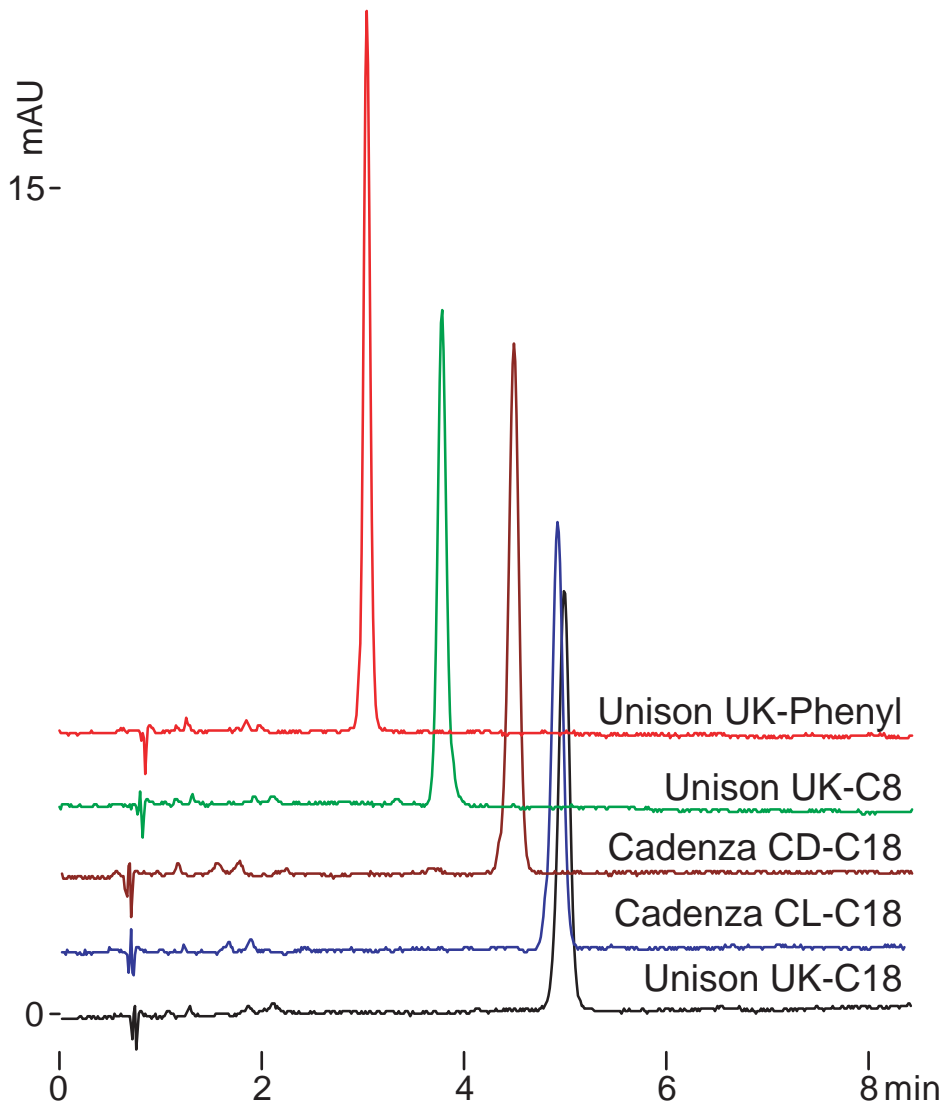

75 x 3 mm, water / acetonitrile / HCOOH = 60 / 40 / 0.1
 0.5 mL/min (5-6MPa), 37 deg.C, 260 nm, 10ppm (5uL)