


Unison UK-Silica

250 x 4.6 mm


Application

Ingredients of Eye-Drops (allantoin and chlorpheniramine maleate)

点眼液成分 (アラントイン, マレイン酸クロルフェニラミン)

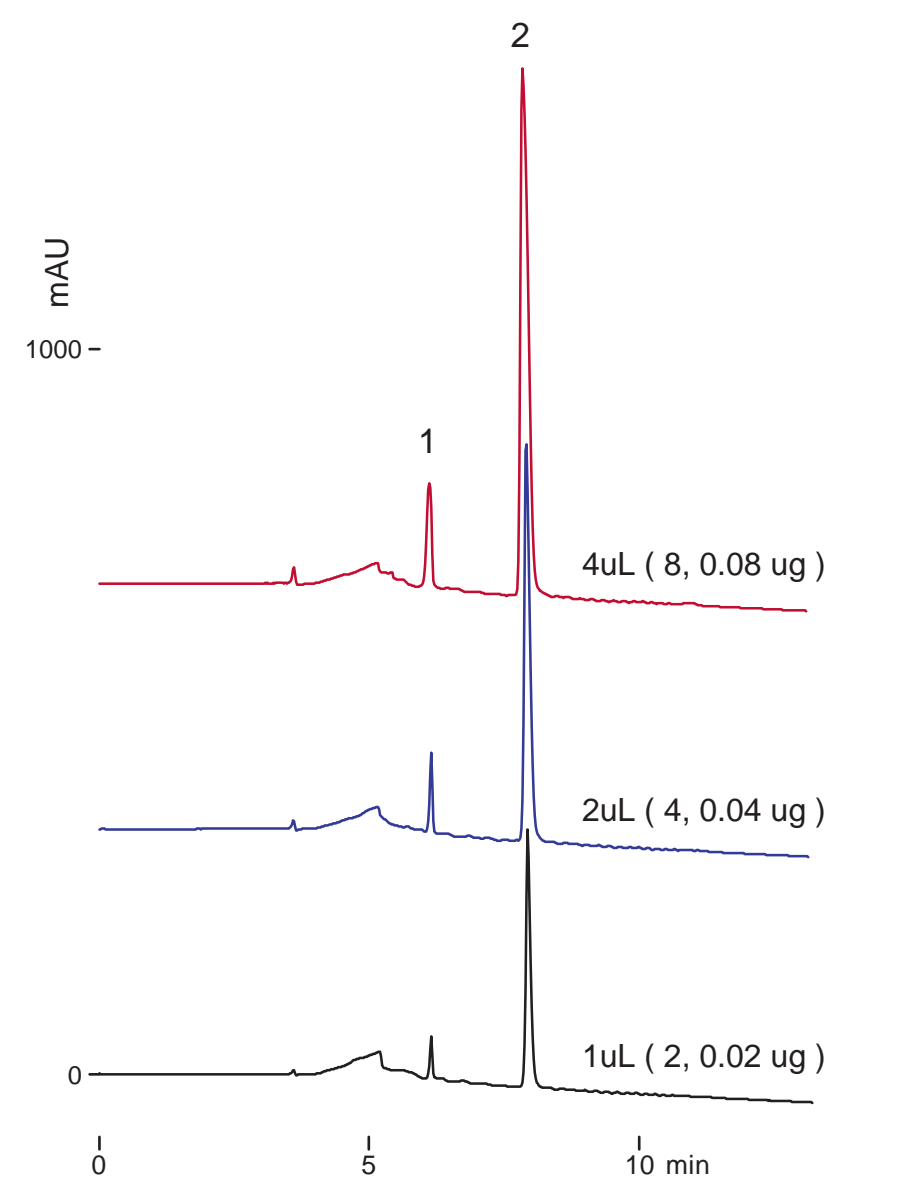


1

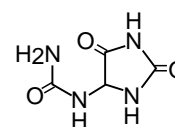

allantoin

2

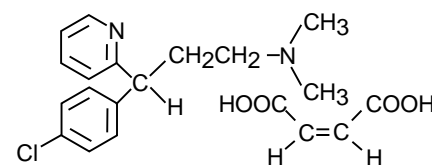

chlorpheniramine, maleate

Unison UK-Silica, 250 x 4.6 mm
 A: acetonitrile / acetic acid = 100 / 0.1
 B: water / acetic acid = 100 / 0.1
 0 - 40%B (0-10min)
 1 mL/min , 37 deg.C, 220 nm