


Unison UK-Amino

250 x 3 mm

Application

### Nucleosides and Bases

核酸塩基・ヌクレオシド

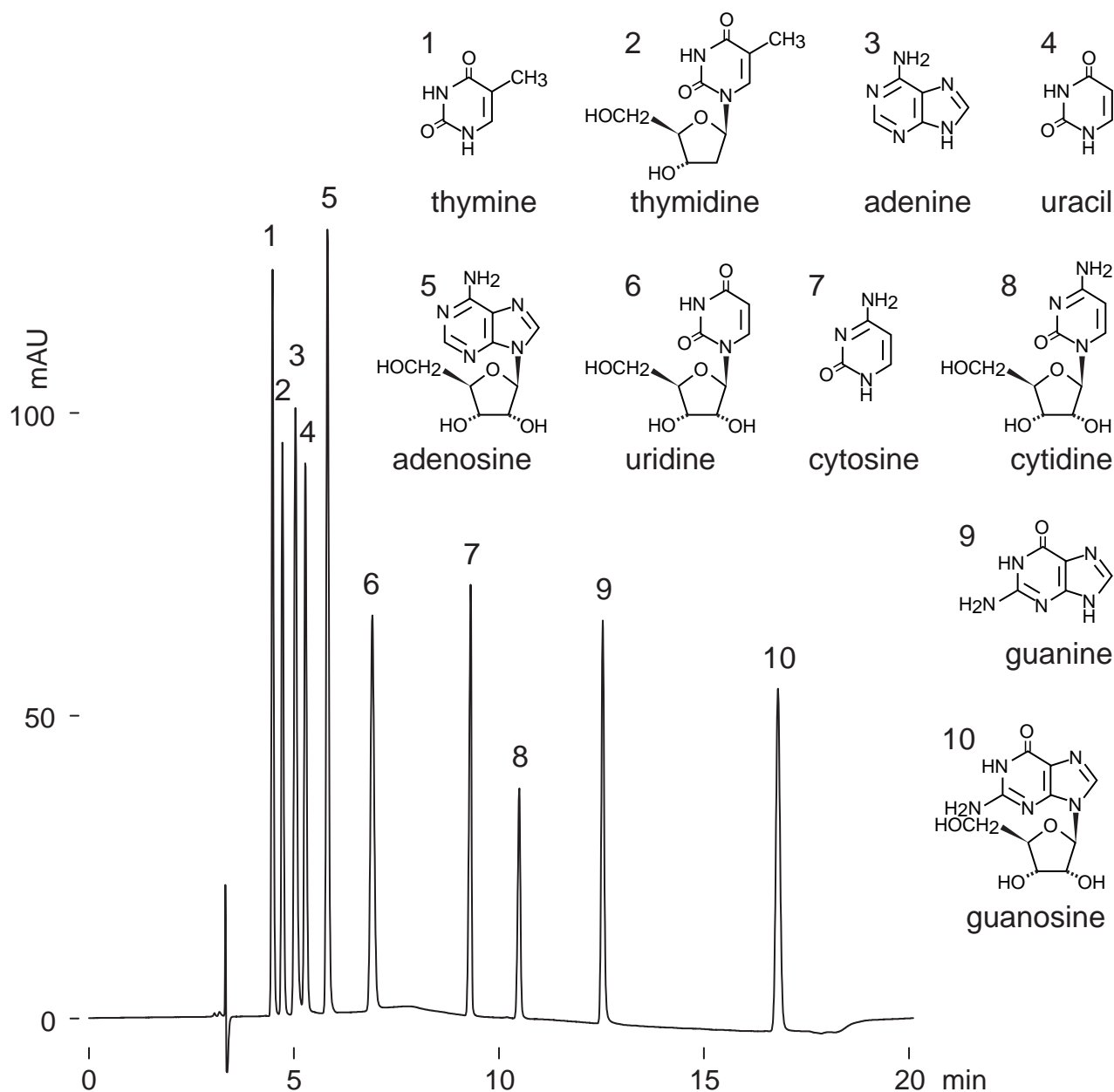

Unison UK-Amino, 250 x 3 mm

A: hexane /tetrahydrofuran /water /trifluoroacetic acid = 5 /91 /4 /0.15

B: tetrahydrofuran /water /trifluoroacetic acid = 50 /50 /0.1

0-33%B(0-15min)

0.4 mL/min (7MPa), 60 deg.C, 260nm

1uL (0.05-0.14ug, diluted with acetonitrile /TFA)