



Unison UK-Phenyl
 Unison UK-C8
 Cadenza CW-C18
 Cadenza CL-C18
 Unison UK-C18
 Cadenza CD-C18

75 x 2 mm


Application

PFAAs


パーフルオロアルキル酸

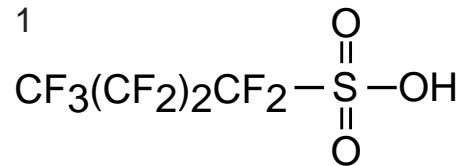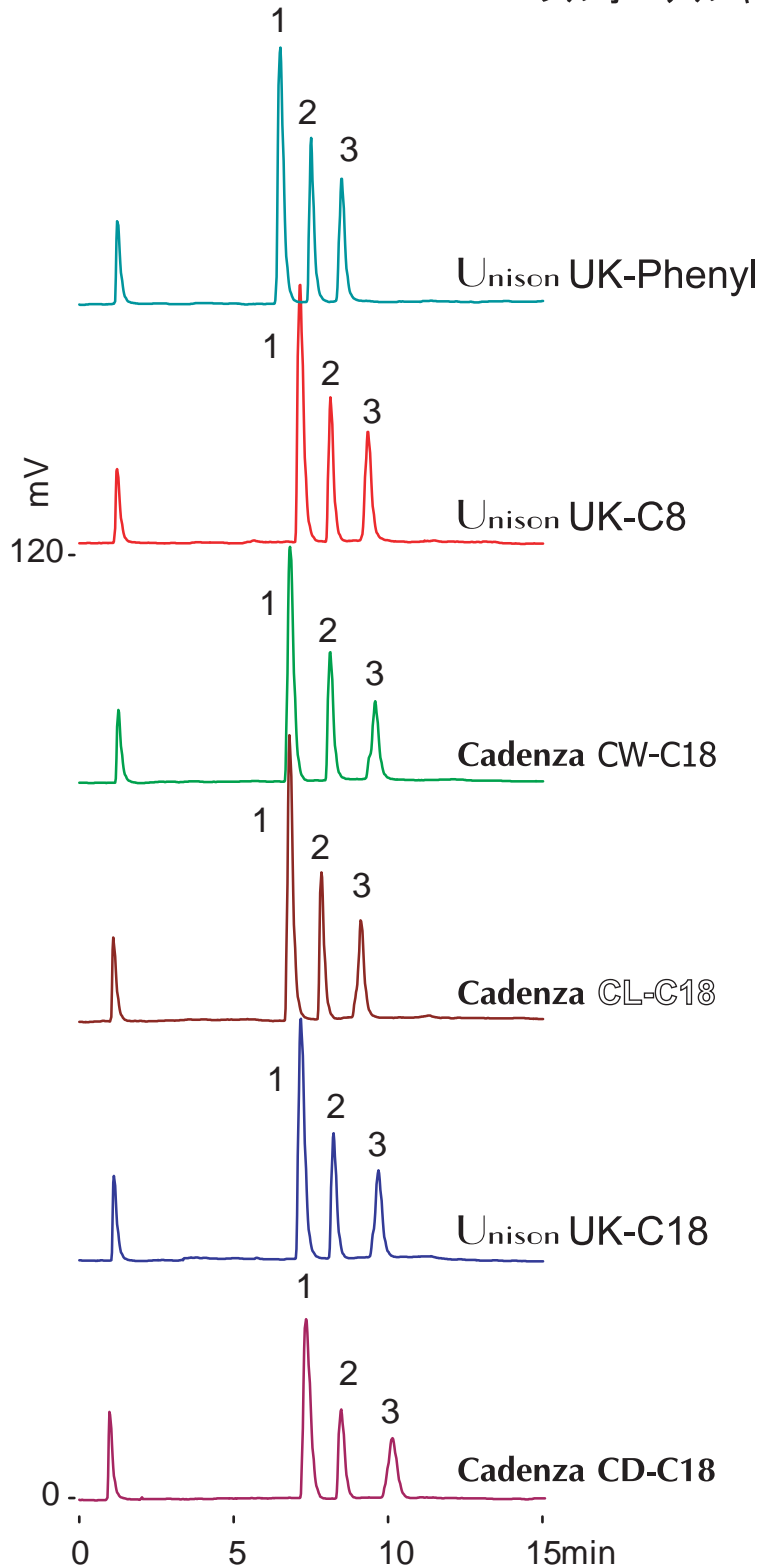

PFBS
 perfluorobutanesulfonic acid

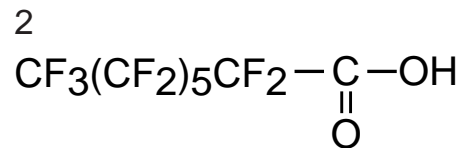

PFOA
 perfluorooctanoic acid

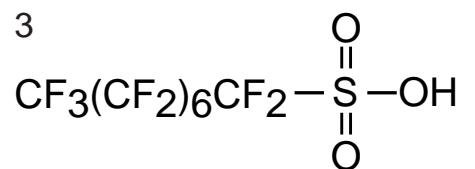

PFOS
 perfluorooctanesulfonic acid

75 x 2 mm,
 A: 10mM ammonium acetate
 B: acetonitrile
 0-70 %B (0-10min)
 0.2 mL/min, 37 deg.C
 ELSD, 3 uL (1.5-1.8ug)