

Unison UK-Amino



250 x 4.6 mm

Application

2-Amino Benzoic Acid (2-AA) Labeled Oligosaccharides

2-アミノ安息香酸(2-AA)ラベル化オリゴ糖

## Fetuin Oligosaccharides

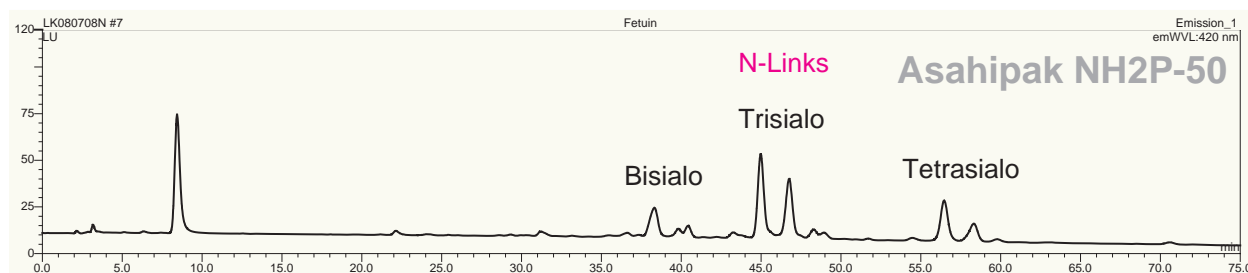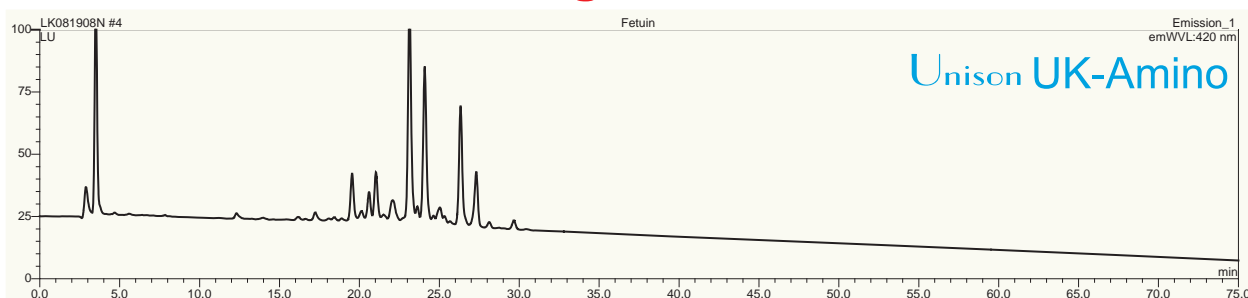

Ref.) Anumula *et.al.*, Glycobiology, vol. 8 no.7 pp.685-694, 1998

Unison UK-Amino, 250 x 4.6 mm  
 A: ACN / AcOH / THF = 97 / 2 / 1  
 B: water / AcOH / TEA / THF = 91 / 5 / 3 / 1  
 30%B (0-2min), 30-95%B (2-80min), 95%B (80-95min)  
 1 mL/min, 50 deg.C. FLD (Ex.360nm, Em.420nm)

Courtesy of Elizabeth Higgins, Glycosolutions Corporation, MA - USA