


Unison UK-Phenyl
 Unison UK-C8
 Unison UK-C18
 Cadenza CD-C18
 Cadenza CW-C18

75 x 3 mm


Application

Malachite Green and Metabolite


マラカイトグリーンと代謝物

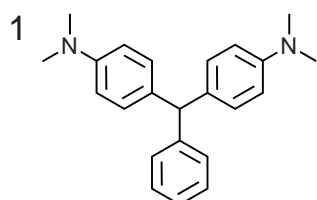

leucomalachite green

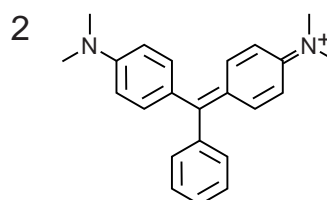

malachite green

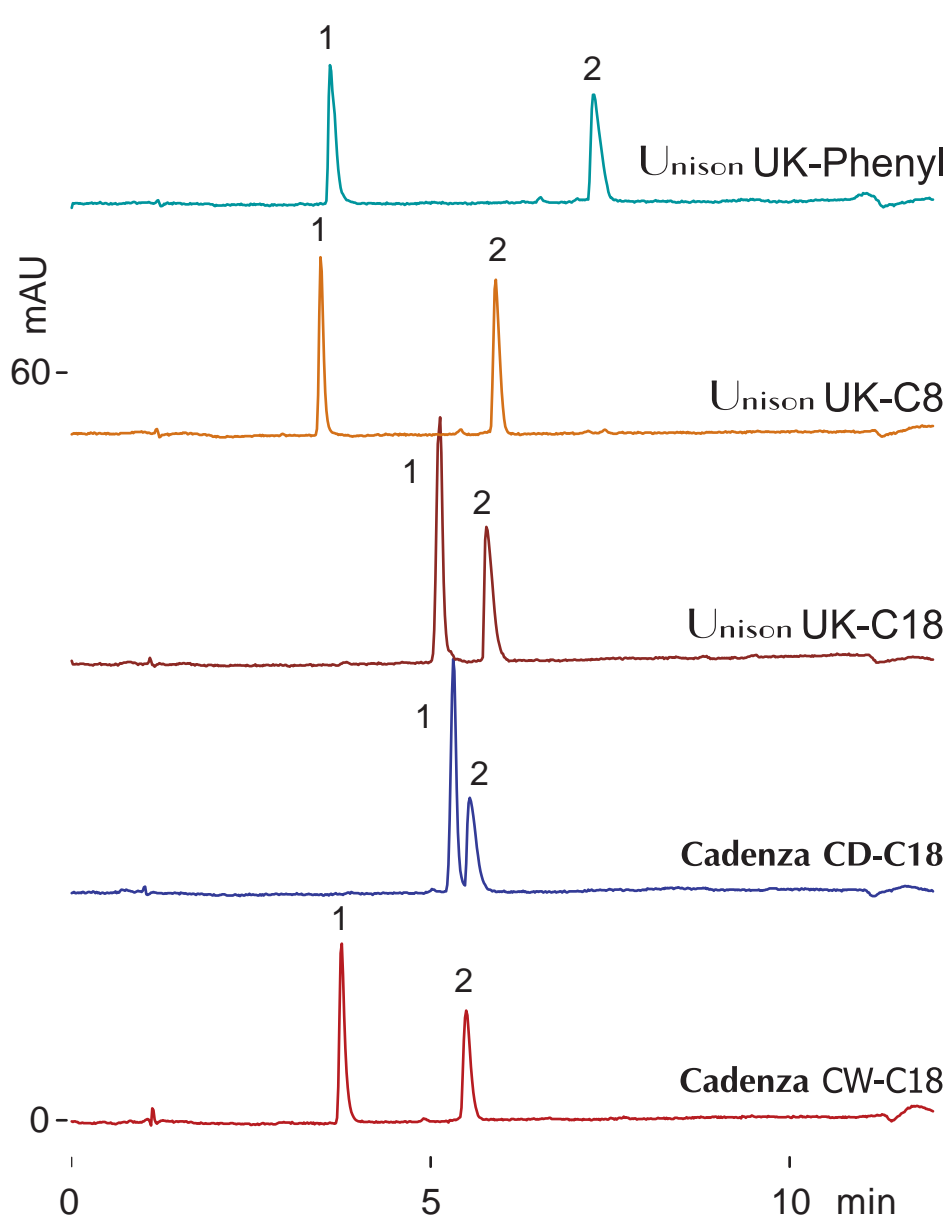

75 x 3 mm
 A: water /formic acid
 = 100 /0.1
 B: acetonitrile /formic acid
 = 100 /0.1
 30 - 70 %B (0-10 min)
 0.4 mL/min (5 MPa)
 37 deg.C, 260 nm
 1 uL (100 ng)