




Unison UK-Amino

50 x 2 mm  
250 x 2 mm


Application

Acidic Amino Acids

酸性アミノ酸

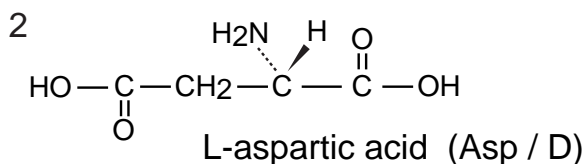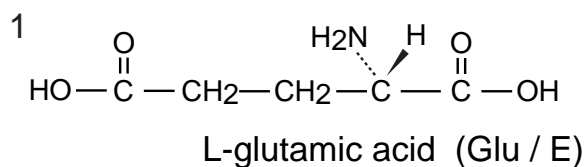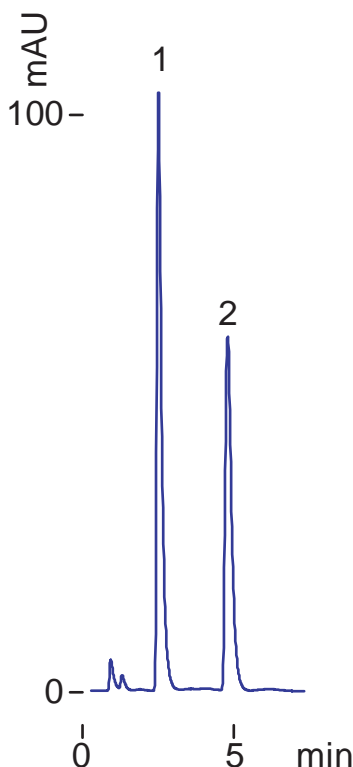

Unison UK-Amino, 50 x 2 mm  
A: acetonitrile  
B: water /acetic acid = 100 /1  
70-90 %B (0-5min)  
0.2 mL/min (3MPa)  
37 deg.C, ELSD, 1 uL (1ug)

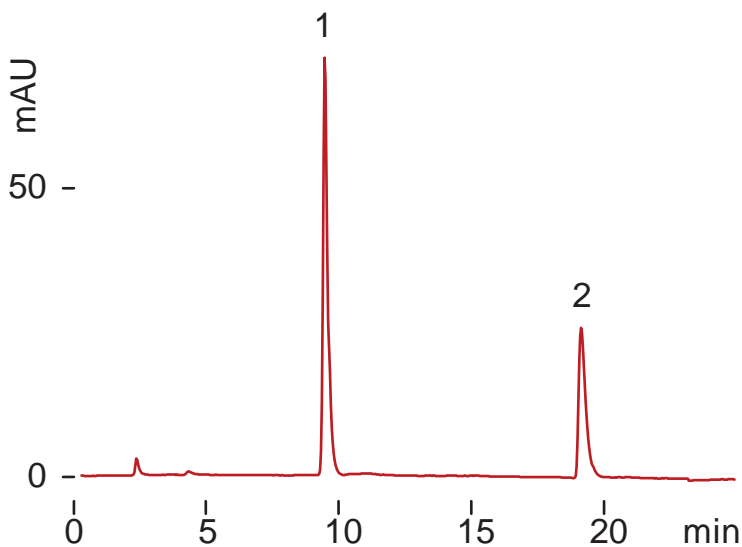

Unison UK-Amino, 250 x 2 mm  
A: acetonitrile  
B: water /acetic acid = 100 /1  
70-90 %B (0-20min)  
0.2 mL/min (16MPa)  
37 deg.C, ELSD, 2 uL (2ug)