
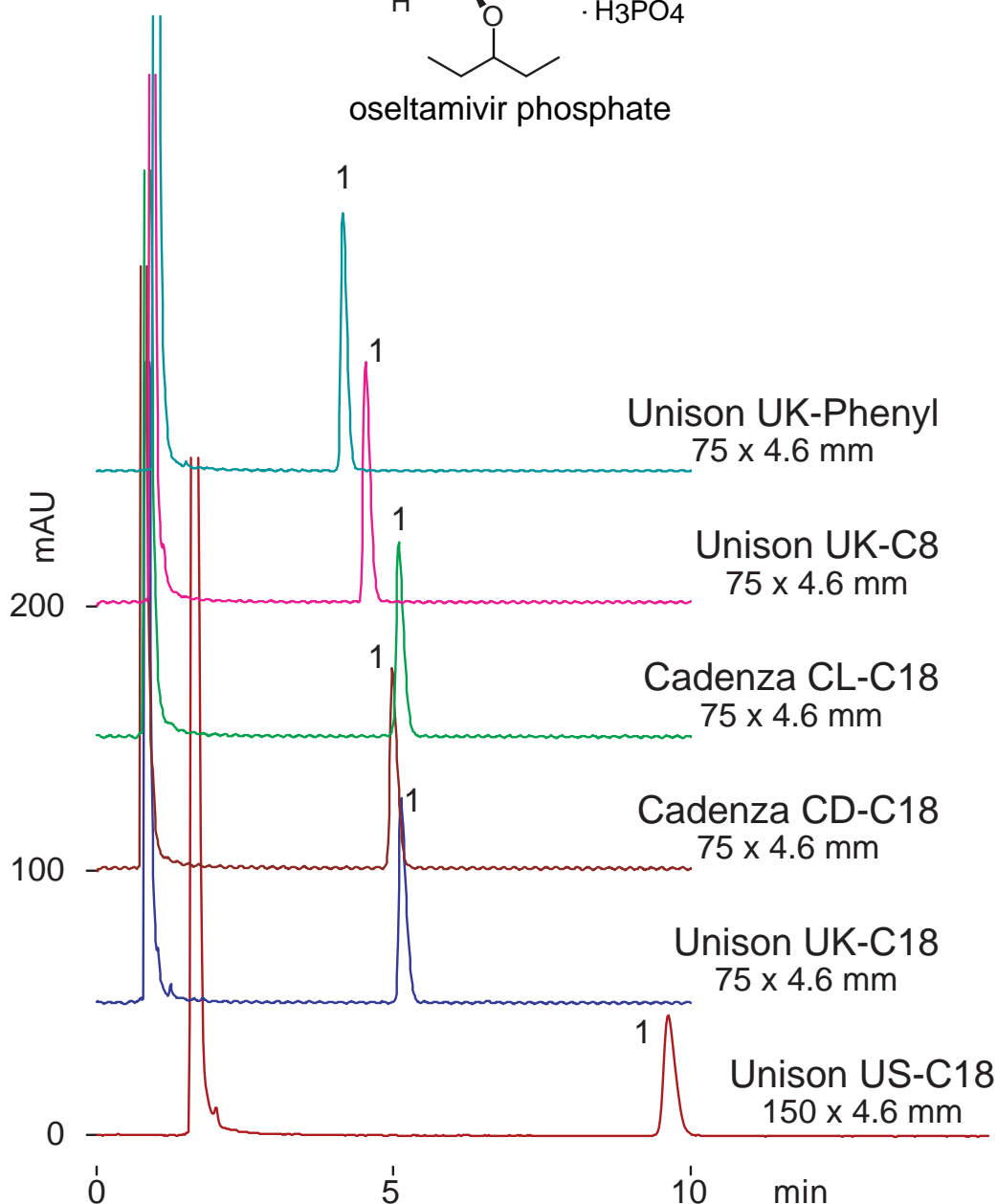


Unison UK-Phenyl Cadenza CD-C18 75 x 4.6 mm
 Unison UK-C8 Unison UK-C18 150 x 4.6 mm
 Cadenza CL-C18 Unison US-C18

Application

Oseltamivir phosphate (Anti-flu medicine)

オセルタミビル リン酸塩 (インフルエンザ治療薬)

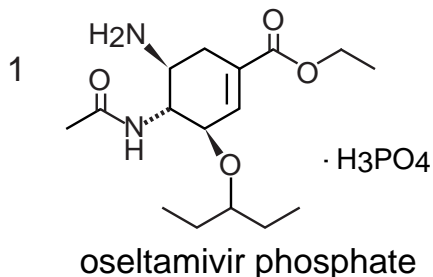

water / acetonitrile / trifluoroacetic acid = 74 / 26 / 0.1
 1 mL/min (4-5MPa), 37 deg.C, 220 nm, 1 uL (1ug, DMSO)