
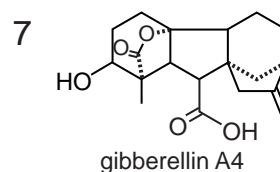
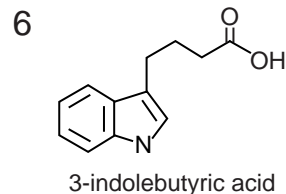
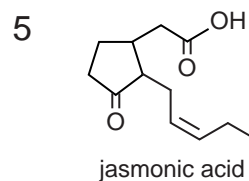
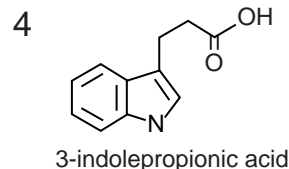
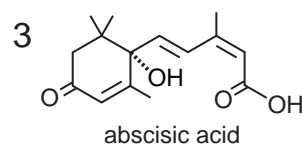
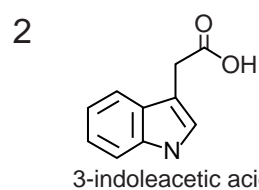
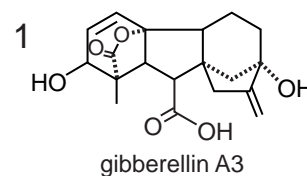


Unison UK-Phenyl
 Unison UK-C8 Cadenza CD-C18
 Cadenza HS-C18 Cadenza CL-C18
 Cadenza CW-C18 Unison UK-C18 150 x 3 mm

Application

Phytohormones

植物ホルモン

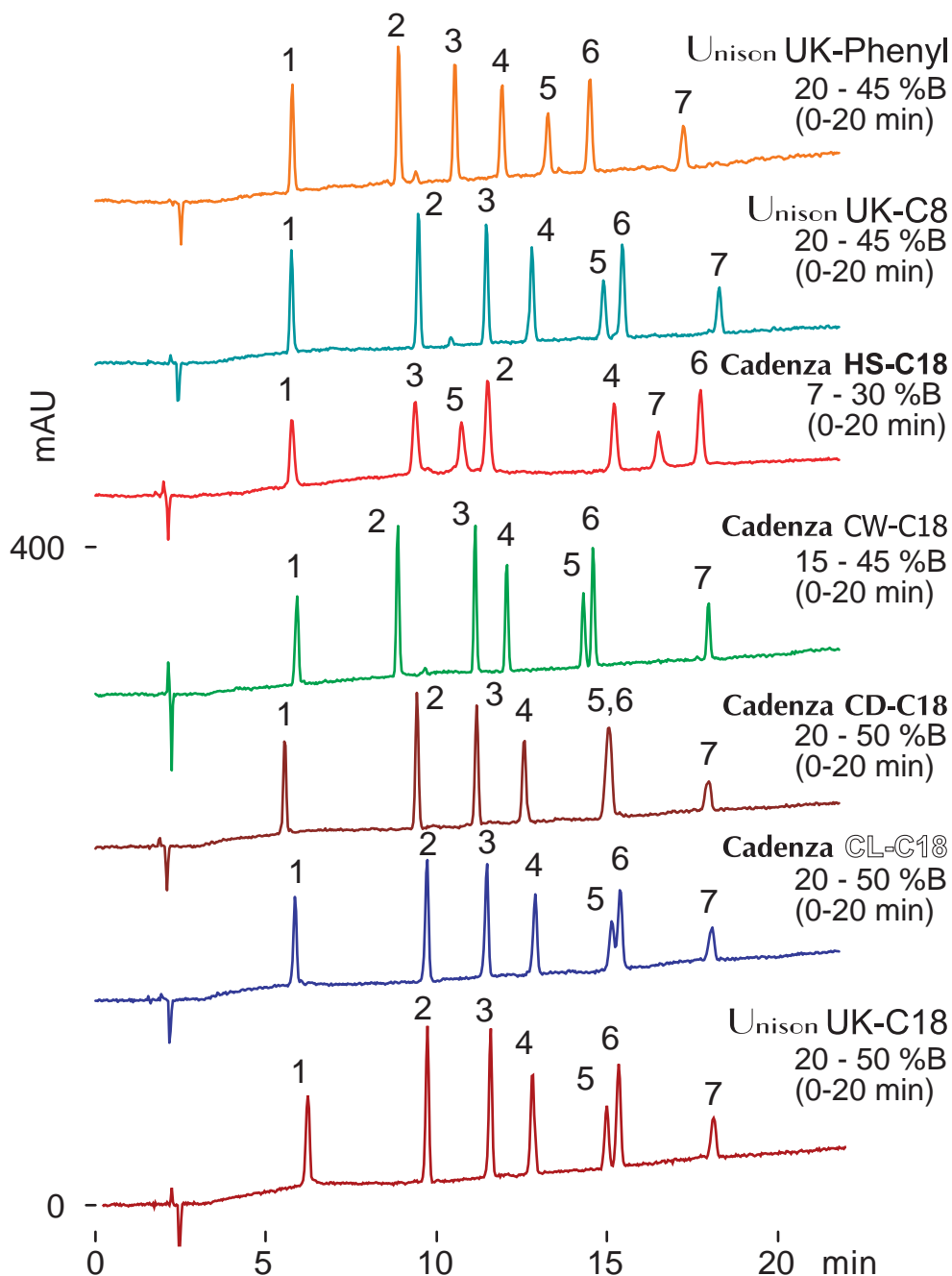

150 x 3 mm

A: water / formic acid = 100 / 0.1

B: acetonitrile / formic acid = 100 / 0.1

0.4 mL/min (8-10 MPa)

37 deg.C, 210 nm, 3 uL (0.3-3 ug)