


**Scherzo SM-C18**  
Unison UK-C18

150 x 3 mm

Application

Acetic acid and trihaloacetic acids

酢酸, トリハロ酢酸

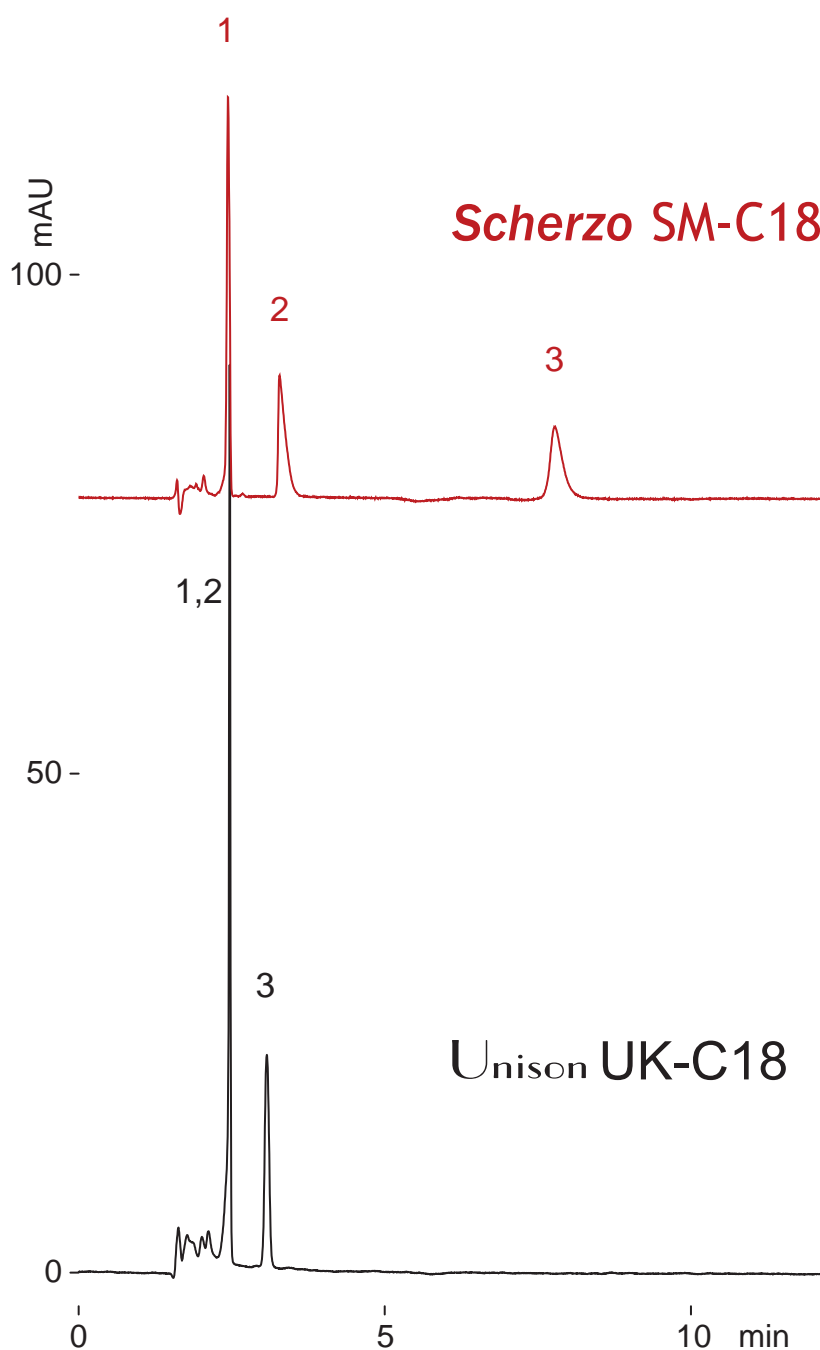

- 1  $\text{CH}_3\text{COOH}$   
acetic acid (100ppm)
- 2  $\text{CF}_3\text{COOH}$   
trifluoroacetic acid (100ppm)
- 3  $\text{CCl}_3\text{COOH}$   
trichloroacetic acid (9ppm)

150 x 3 mm  
100mM  $\text{NH}_4\text{H}_2\text{PO}_4$  / acetonitrile = 80 / 20  
0.4 mL/min (9-10MPa), 37 deg.C, 210 nm  
20 uL