


Presto FF-C18
Cadenza CD-C18

250 x 4.6 mm


Application

Diisononyl phthalate


フタル酸ジイソノニル

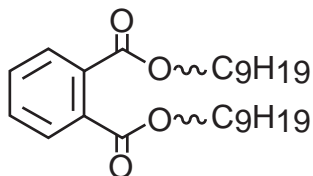

diisononyl phthalate
DINP

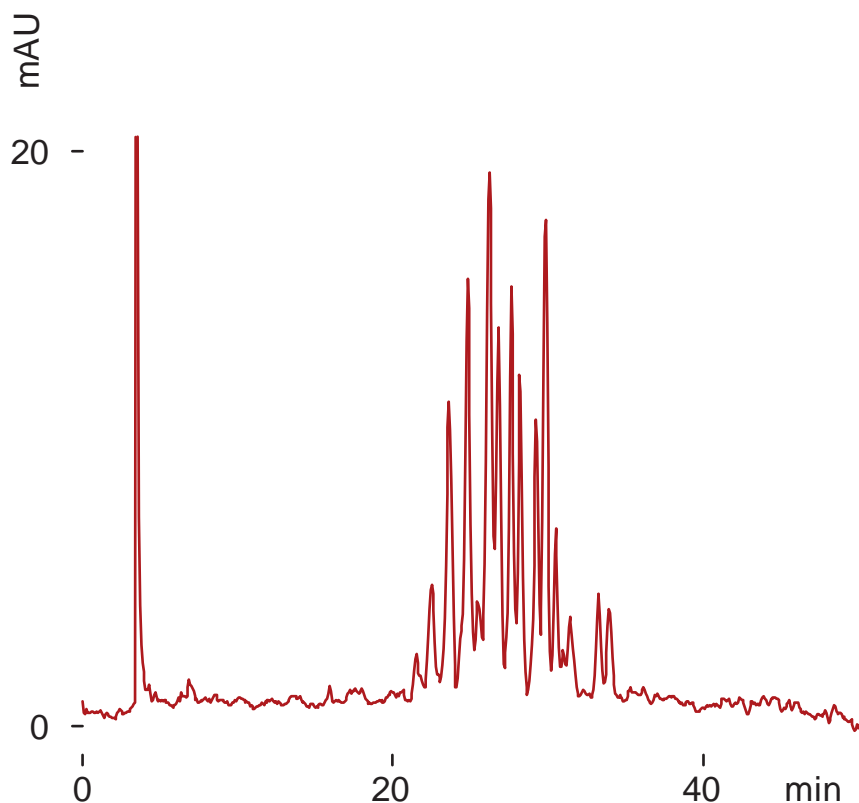

Presto FF-C18
250 x 4.6 mm
water / acetonitrile = 40 / 60
0.45 mL/min (28MPa)
37 deg.C, 230 nm, 3 uL (3ug)

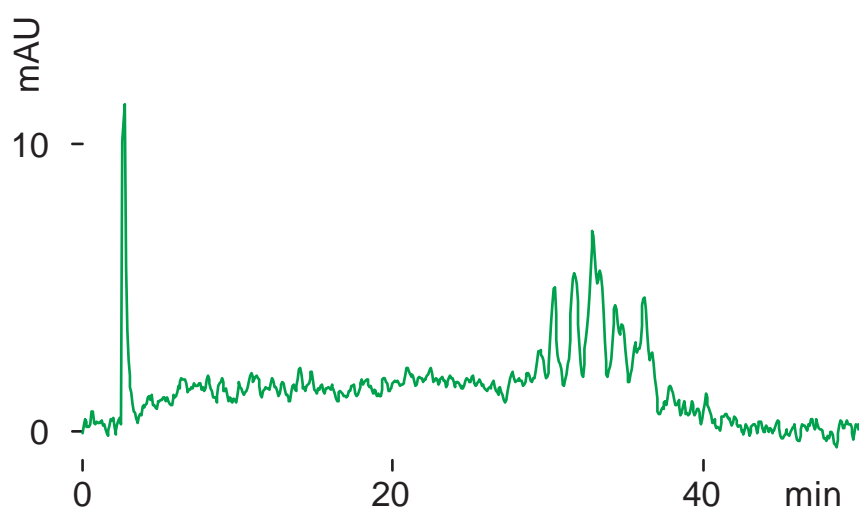

Cadenza CD-C18
250 x 4.6 mm
water / acetonitrile = 10 / 90
1 mL/min (9MPa)
37 deg.C, 230 nm, 3 uL (3ug)