

Cadenza HS-C18  
Cadenza CW-C18  
Unison UK-Phenyl  
Unison UK-C8



Scherzo SM-C18  
Cadenza CL-C18  
Cadenza CD-C18  
Unison UK-C18

75 x 4.6 mm

Application

### Bisphenol A

ビスフェノールA

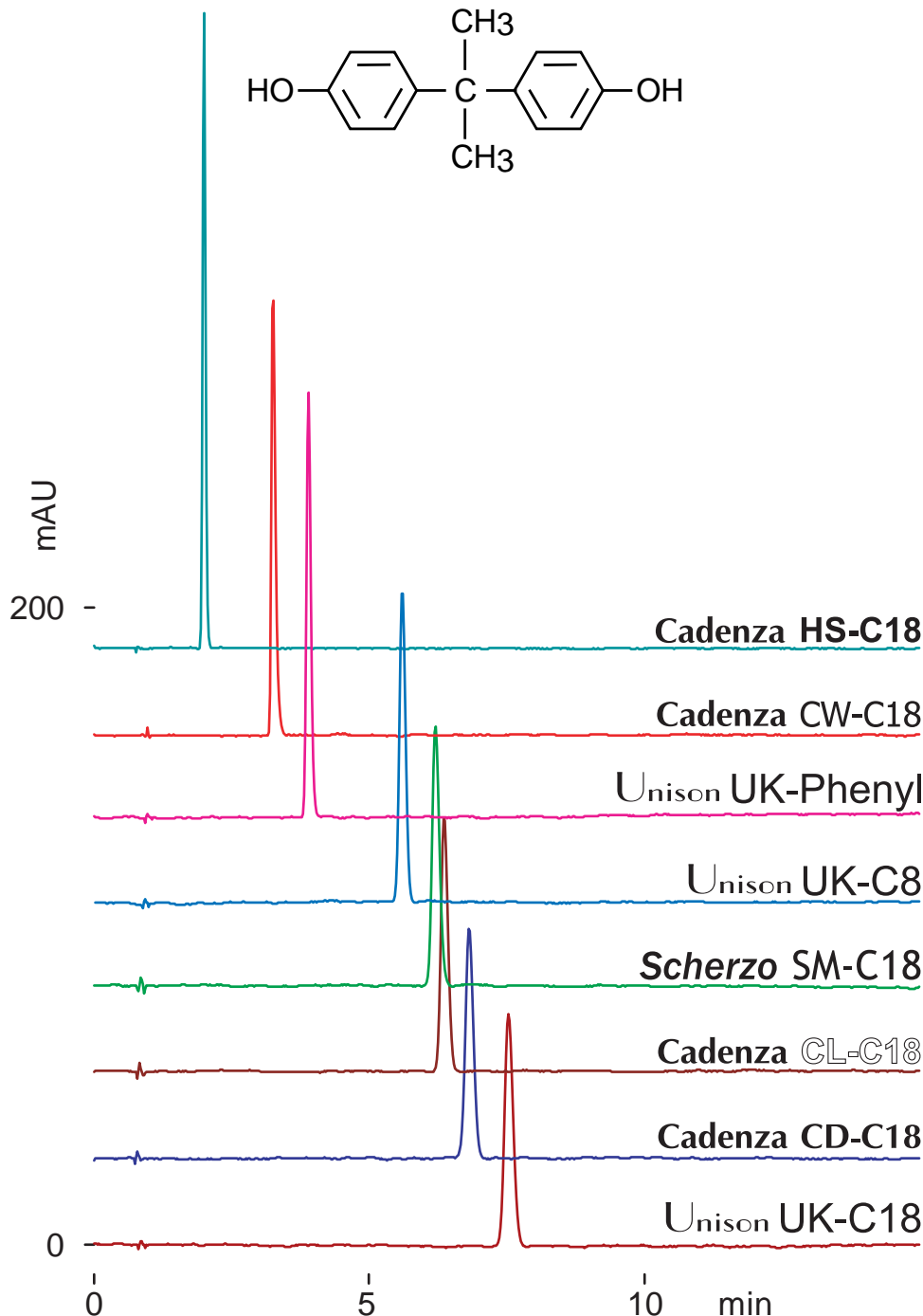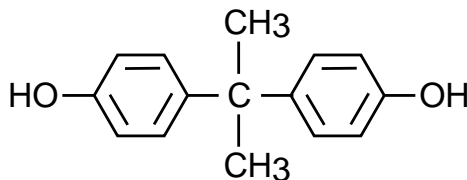

75 x 4.6 mm  
water / acetonitrile / acetic acid = 65 / 35 / 0.1  
1 mL/min (6-7MPa), 37 deg.C, 260 nm, 1 uL (5ug)