

Cadenza HS-C18
Cadenza CW-C18
Unison UK-Phenyl
Unison UK-C8



Cadenza CD-C18
Cadenza CL-C18
Unison UK-C18
Scherzo SM-C18

150 x 4.6 mm

Application

Cetirizine (Antihistamine)

セチリジン(抗ヒスタミン薬)

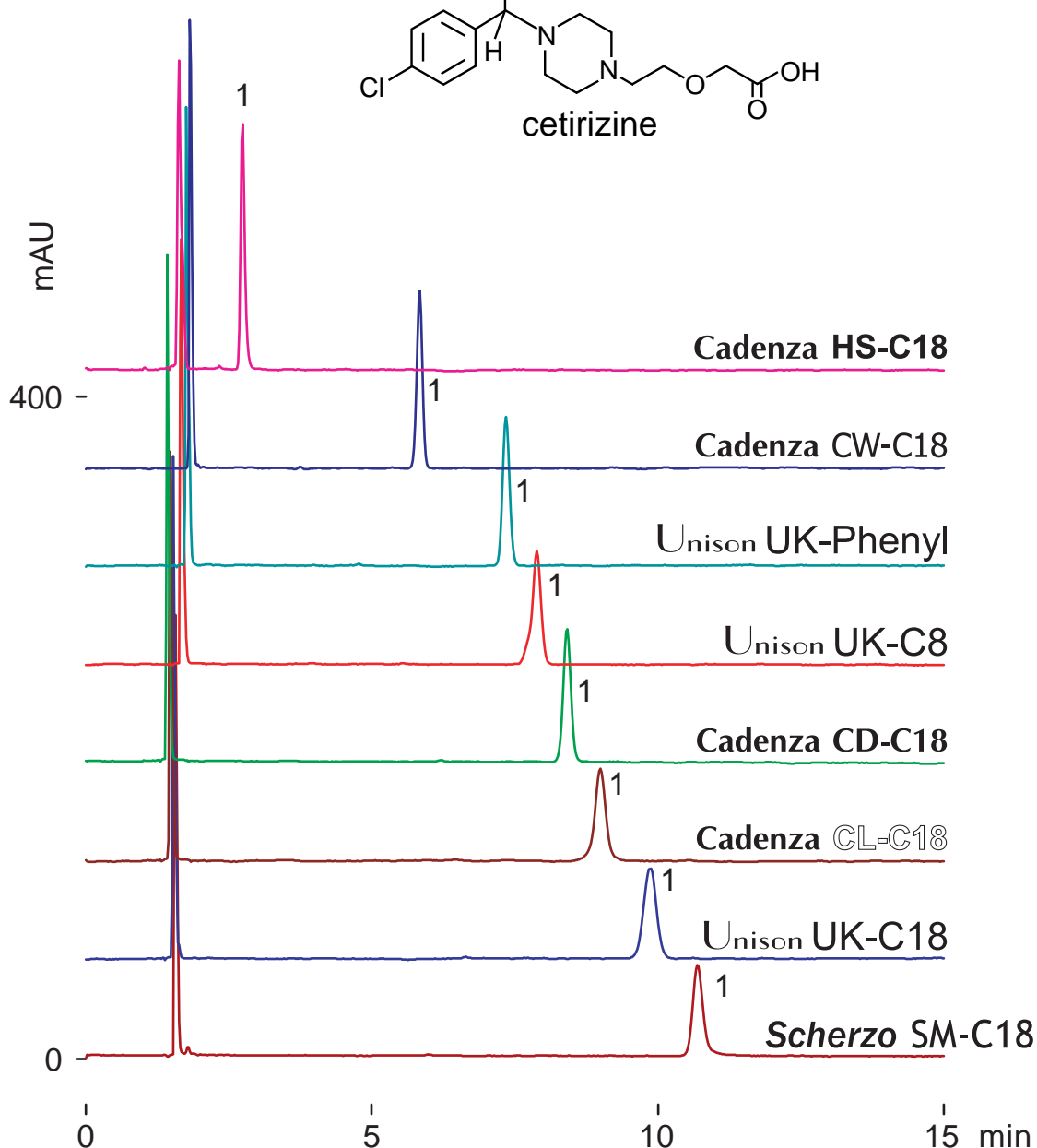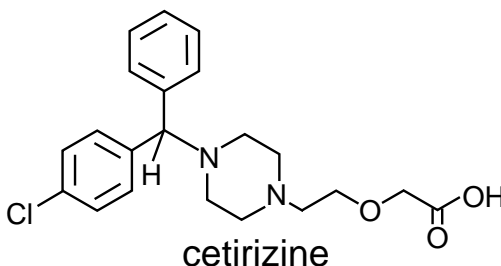

150 x 4.6 mm

100mM ammonium formate / acetonitrile = 67 / 33

1 mL/min (9-12MPa), 37 deg.C, 240 nm

2 uL (2 ug, DMSO)