


Scherzo SM-C18

150 x 2 mm

Application

## Analysis of Water-soluble Vitamins using LCMS-2020

LCMS-2020を用いた水溶性ビタミンの分析

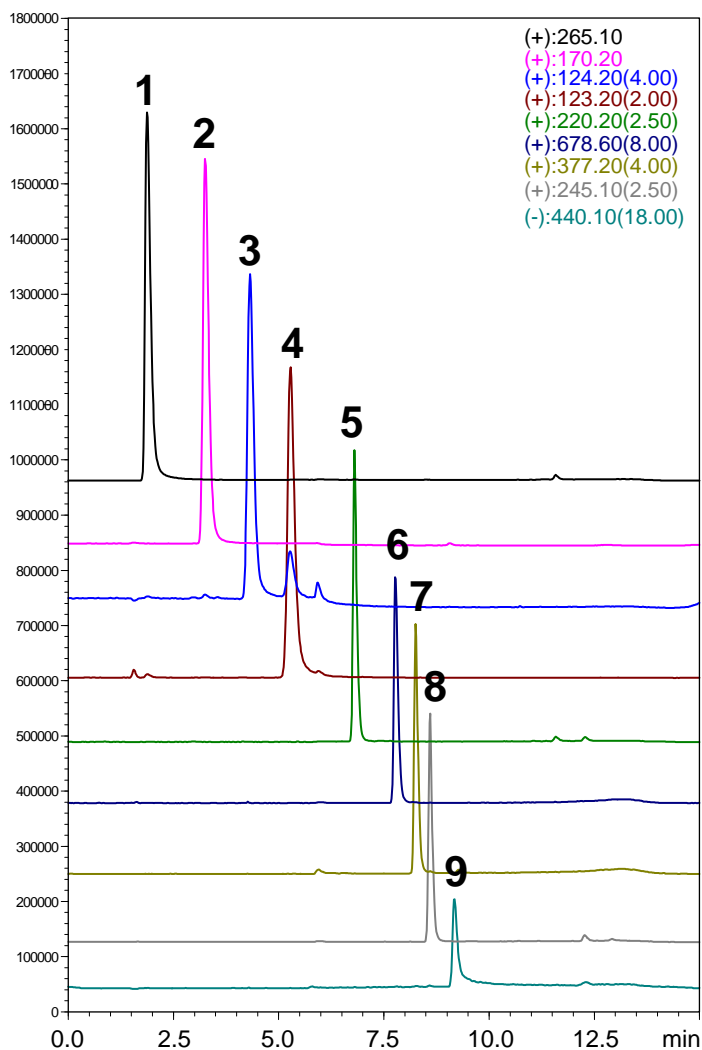

**Shimadzu LCMS-2020**  
(ESI, positive/negative, SIM)

Scherzo SM-C18, 150 x 2 mm  
 A: 5 mmol/L ammonium formate  
 + 0.1% formic acid - water  
 B: acetonitrile  
 0 - 55 %B (0-10 min), 0 %B(10.01-20 min)  
 1,000 pg/uL, 1uL

