






Cadenza CD-C18

150 x 2 mm

Application

LC-MS/MS analysis of nitrofurans metabolites (AOZ, AMOZ, AHD)

ニトロフラン類代謝物 (AOZ, AMOZ, AHD) のLC-MS/MS分析(告示試験法)

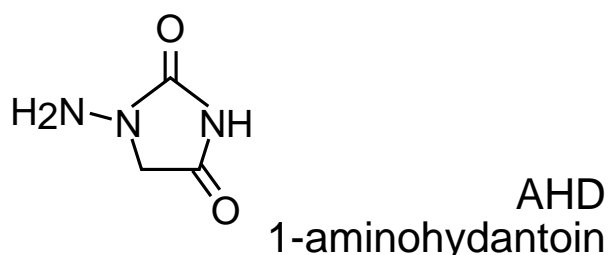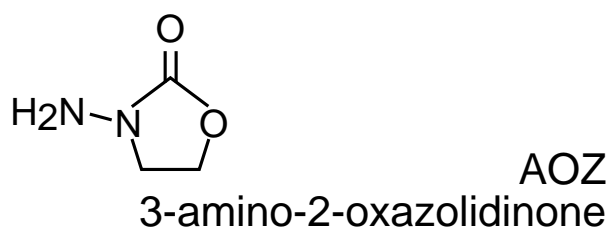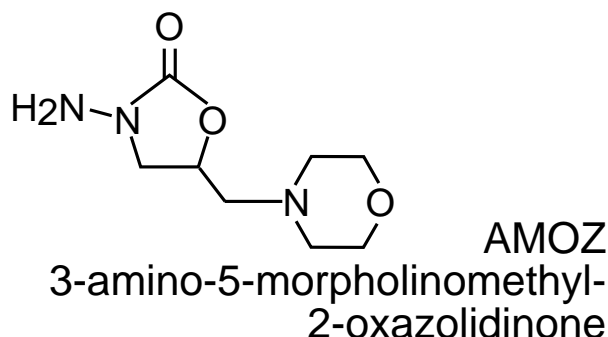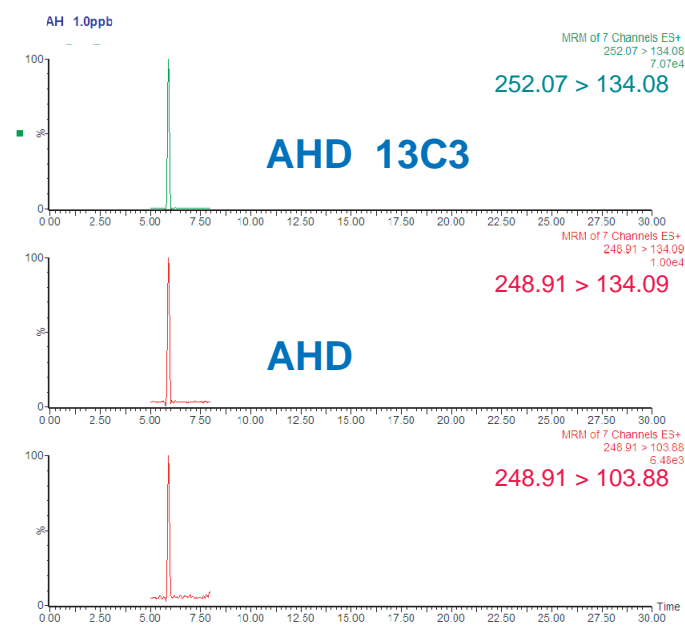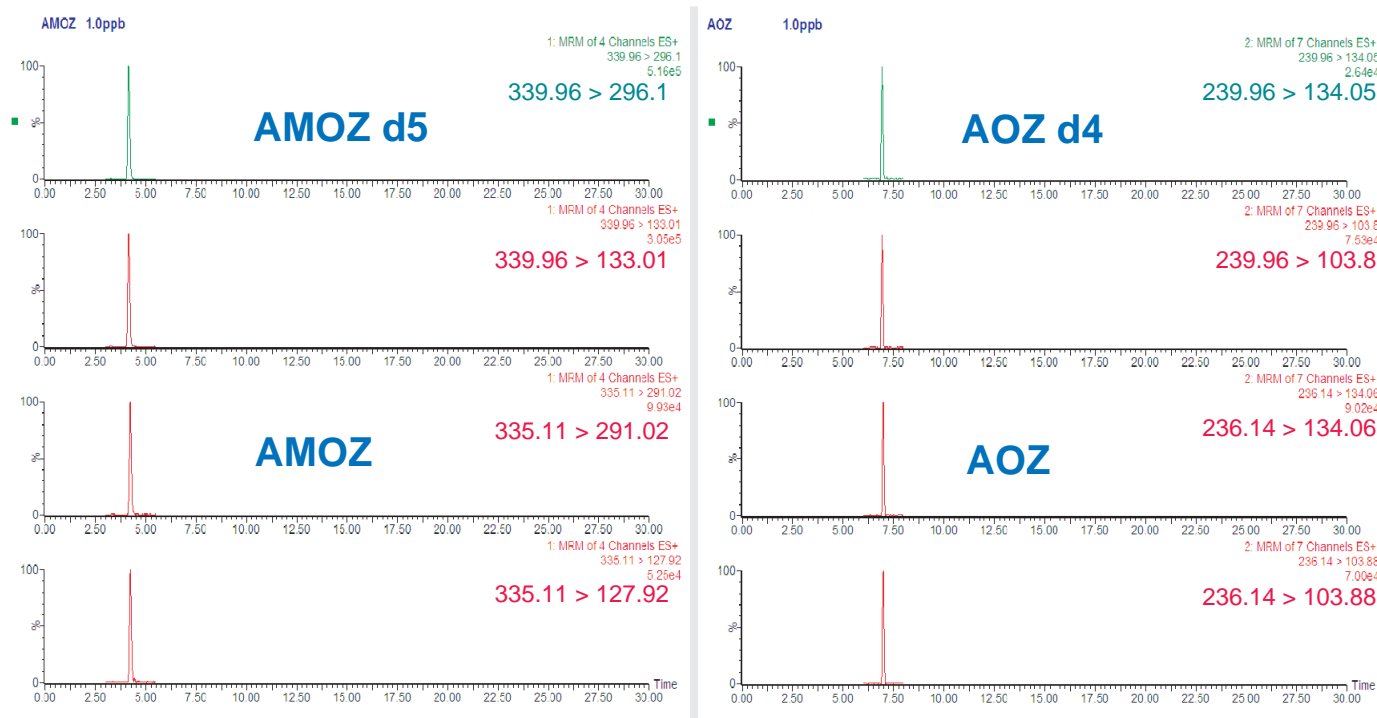

Cadenza CD-C18, 150 x 2 mm
A: 0.1% AcOH
B: acetonitrile
20-80%B (0-15min), Runtime 19min
0.2 mL/min, 40°C
5uL (AOZ, AMOZ, AHD 1.0 ppb, IS 5.0ppb)
ESI, MRM, positive

Courtesy of MAEDA Naoyuki, ASSOCIATION OF MEAT SCIENCE & TECHNOLOGY INSTITUTE, JAPAN