



Presto FF-C18
Cadenza CD-C18

150 x 3 mm

Application

Peptide mapping using LC-FTMS

LC-FTMSによるペプチドマッピング

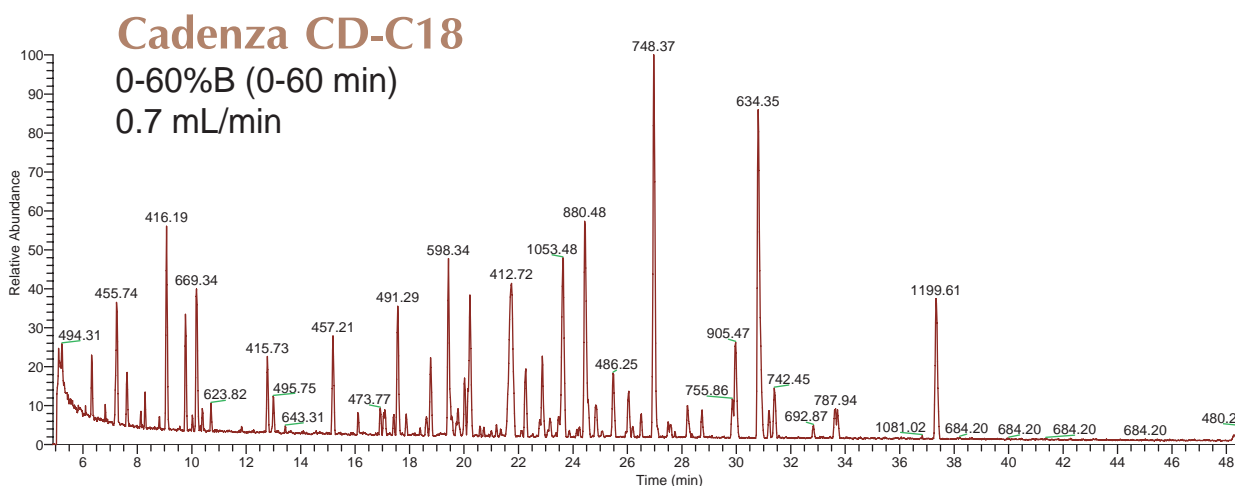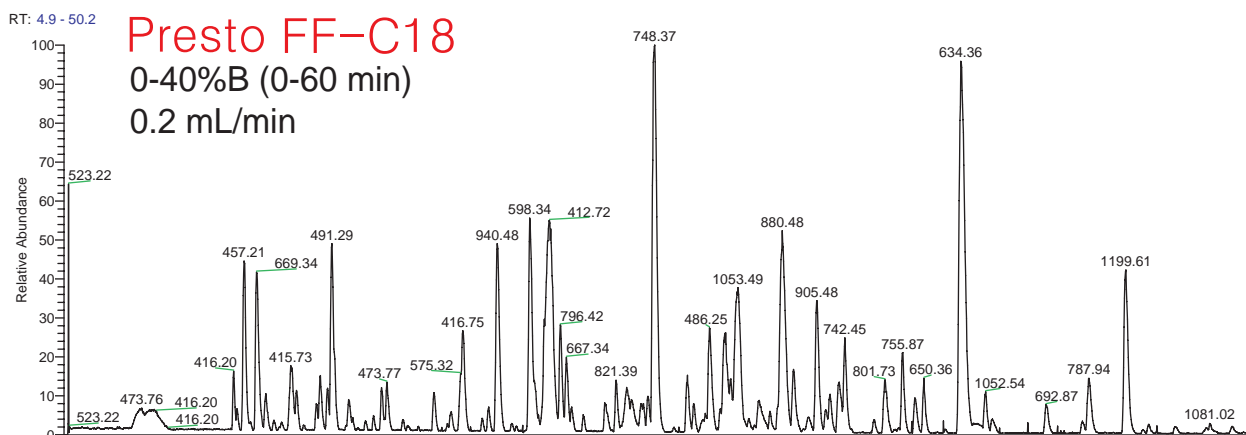

Presto FF-C18, Cadenza CD-C18, 150 x 3 mm

A: water / formic acid / TFA = 100 / 0.5 / 0.02

B: acetonitrile / formic acid / TFA = 100 / 0.5 / 0.02

50 deg.C, 2uL (tryptic digest of alpha-casein)

FTMS (Thermo Scientific Exactive), ESI-positive (Full ms)

Courtesy of Dr. KUWATA Keiko, Kyoto University, JAPAN