


Unison UK-C18
 Scherzo SW-C18
 Scherzo SM-C18
 Scherzo SS-C18

250 x 3 mm

Application

Sialylglycopeptide in egg yolk extract

卵黄抽出物中のシアリルグリコペプチド

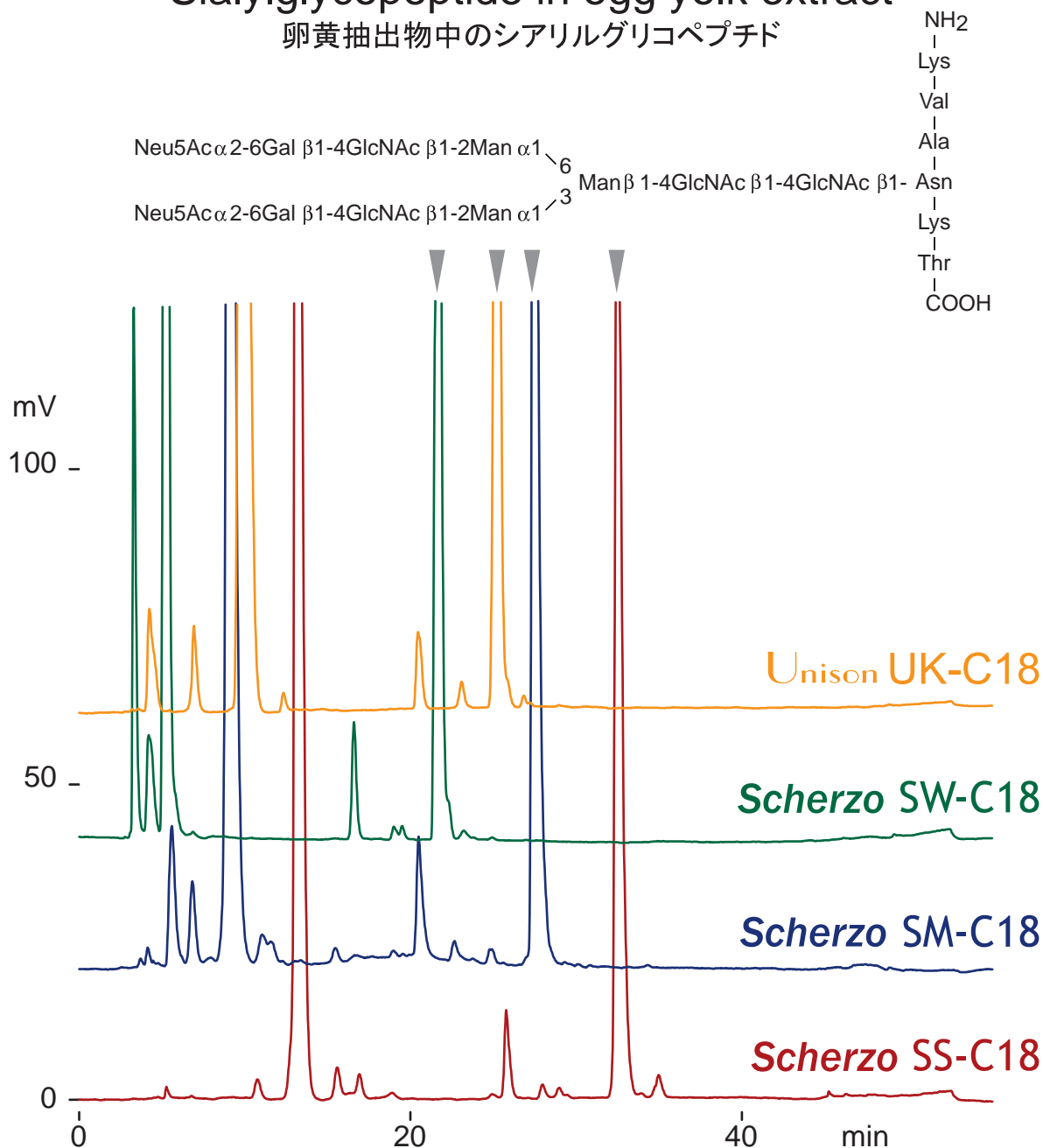

250 x 3 mm
 A: 20mM ammonium acetate
 B: 100mM ammonium acetate /acetonitrile = 50 /50
 0-10 %B(0-40min), 10-100 %B(40-50min)
 0.4 mL/min (12MPa), 37 deg.C
 ELSD (spray chamber 40 deg.C, drift tube 100 deg.C)
 5 uL

Courtesy of FUSHIMI Pharmaceutical Co., LTD.