



Scherzo SS-C18
 Scherzo SM-C18
 Scherzo SW-C18
 Unison UK-C18

150 x 3 mm

Application

Arsenobetaine

アルセノベタイン

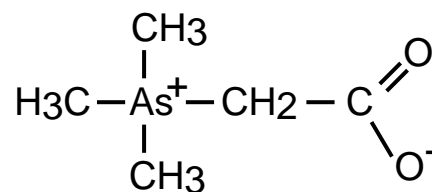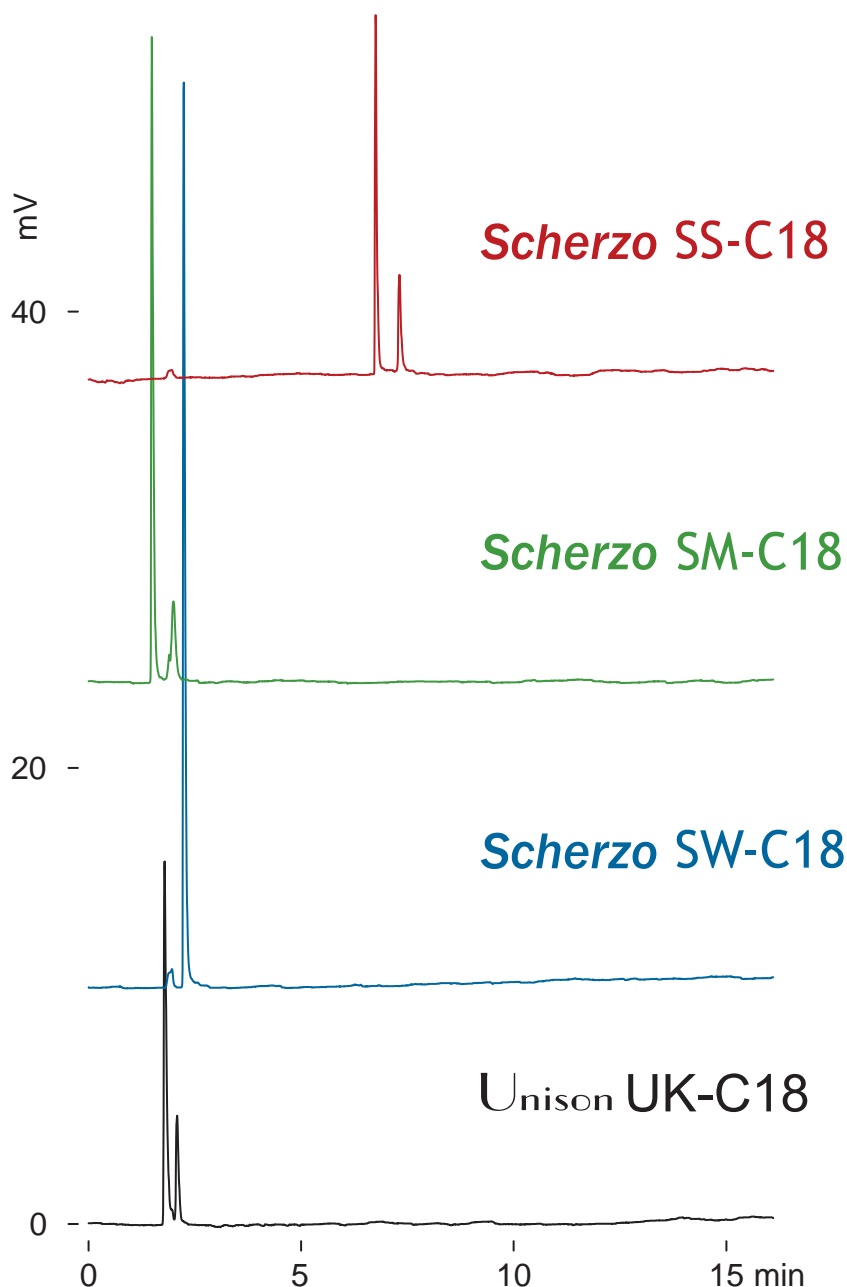

arsenobetaine

150 x 3 mm

A: water / formic acid = 100 / 0.2

B: 100 mM ammonium formate / acetonitrile = 50 / 50

5-80 %B (0-15 min)

0.4 mL/min (8-9 MPa), 37 deg.C, ELSD (spray chamber 50 deg.C, drift tube 100 deg.C)

20 uL (0.49 ug)