






Scherzo SS-C18
 Scherzo SM-C18
 Scherzo SW-C18
 Unison UK-C18

150 x 3 mm

Application

Shikimic acid pathway

シキミ酸経路

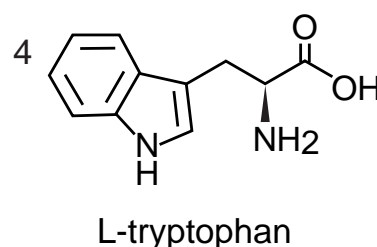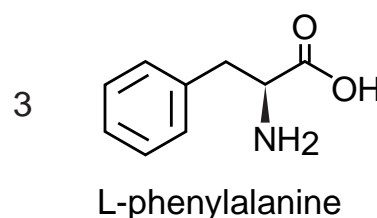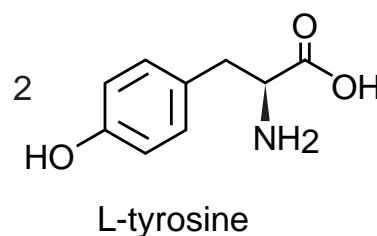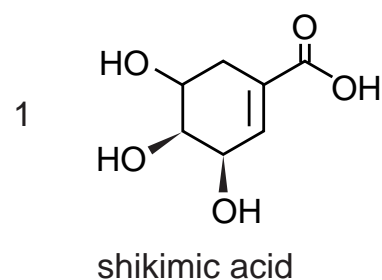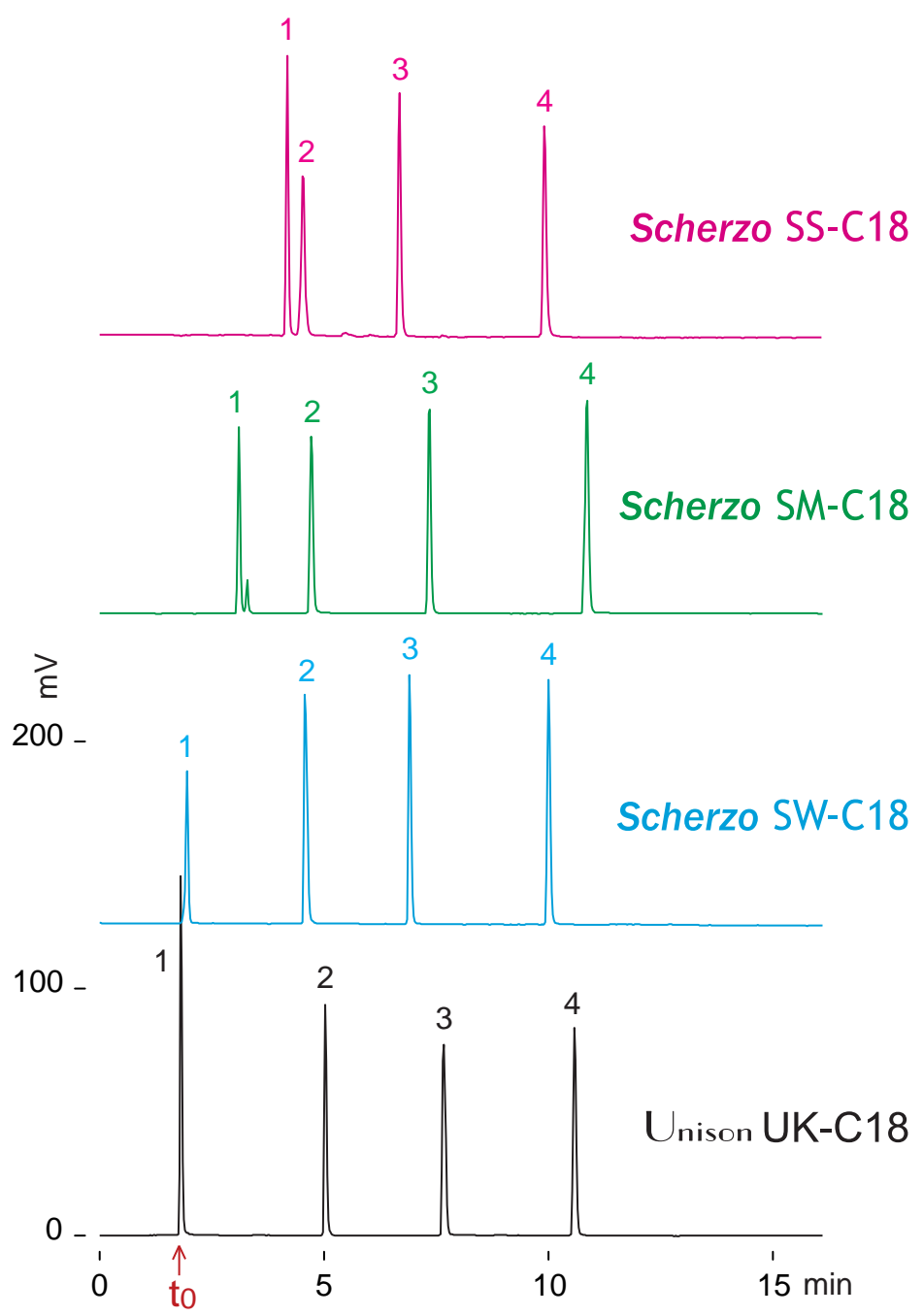

150 x 3 mm

A: 5 mM ammonium acetate

B: 50 mM ammonium acetate / acetonitrile = 90 / 10

0-100 %B (0-10 min)

0.4 mL/min (11 MPa), 37 deg.C, ELSD (spray chamber 50 deg.C, drift tube 100 deg.C)

1 uL (2 ug, 0.06 N HCl aq.)