


**Intrada Amino Acid**

150 x 3 mm

Application

## Batch-to-Batch Reproducibility (Intrada Amino Acid)

バッチ間再現性 (Intrada Amino Acid)

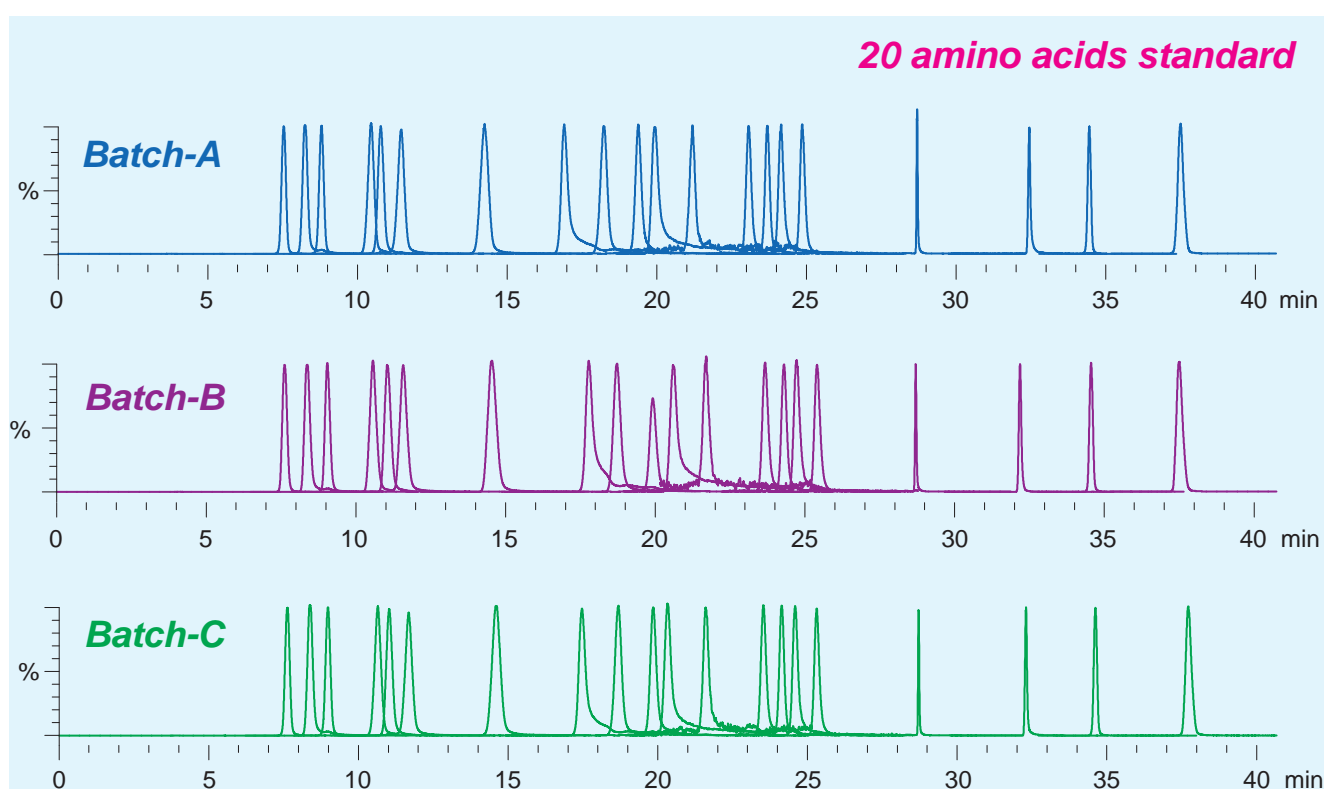

**Intrada Amino Acid, 150 x 3 mm**

**A: ACN / THF / 25mM HCOONH<sub>4</sub> / HCOOH = 9 / 75 / 16 / 0.3 (v/v/v/v)**

**B: ACN / 100mM HCOONH<sub>4</sub> = 20 / 80 (v/v)**

**0 %B (0-12 min), 0-17 %B (12-26 min), 100 %B (26-45 min)**

**0.4 mL/min (7MPa), 35 deg.C, 1.2 uL (1 umol/mL in 0.1N HCl)**

**ESI, positive**