







Scherzo SM-C18

100 x 3 mm

Application

LC-MS: Bisphosphonates (bone resorption inhibitors)

ビスホスホネート系骨粗鬆症治療薬

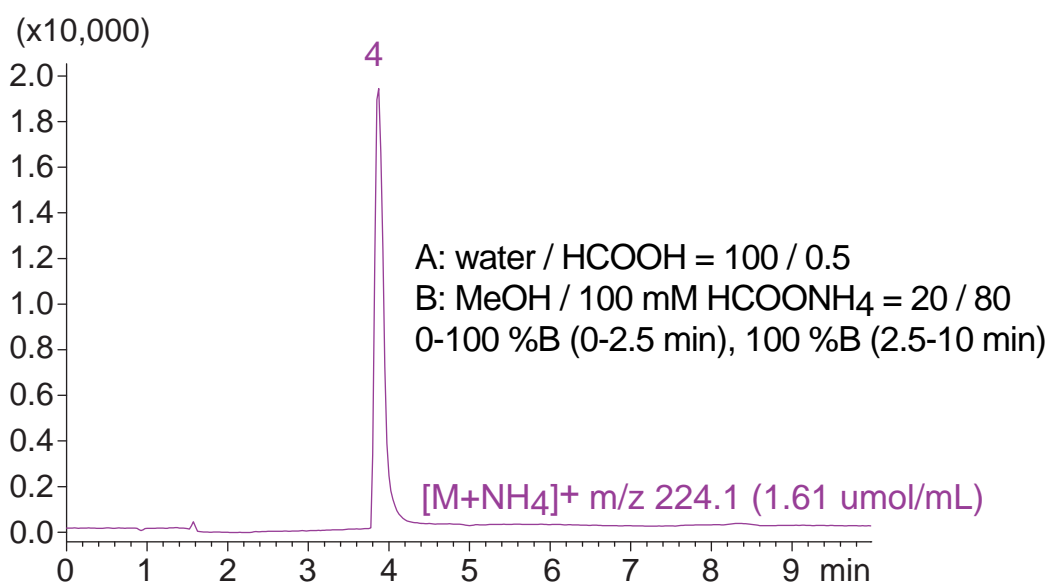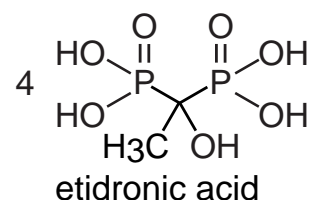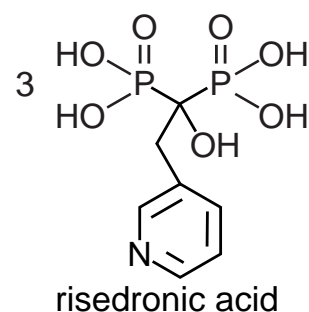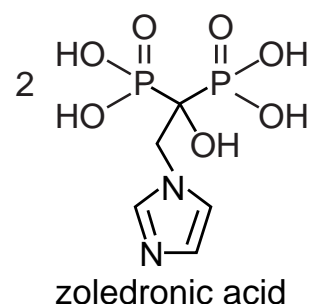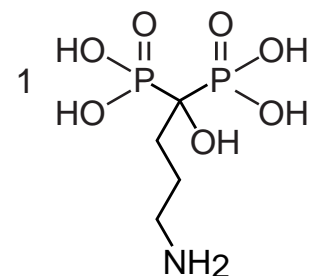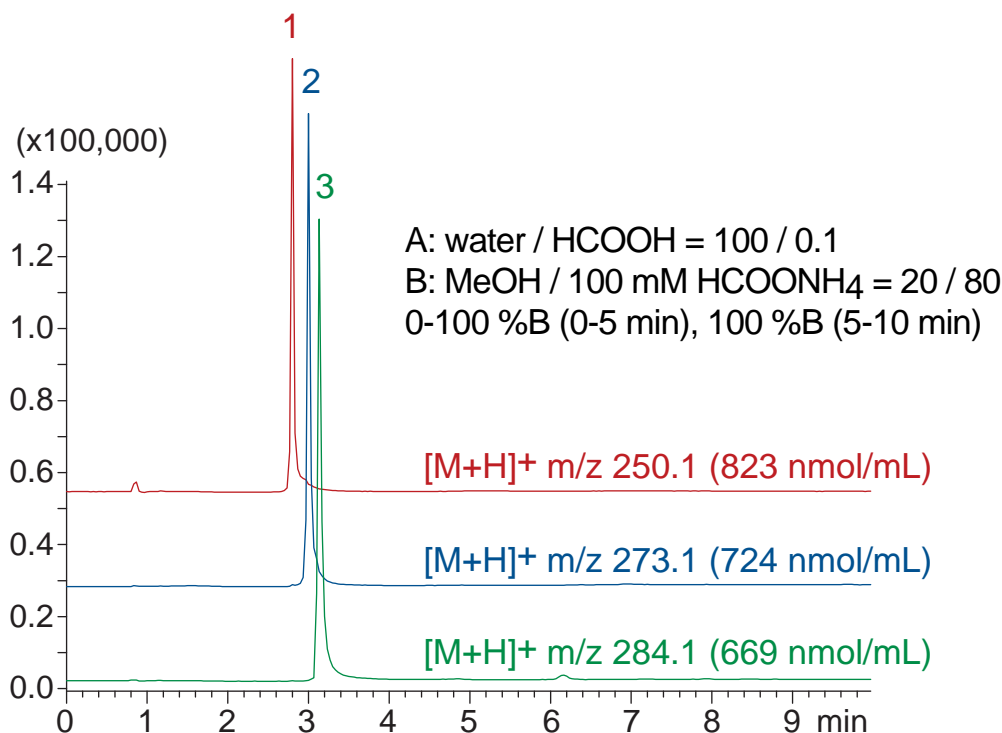

Scherzo SM-C18, 100 x 3 mm
 0.4 mL/mL (8 MPa), 37 °C, 5 uL

INSTRUMENT:

Shimadzu Nexera UHPLC System - Shimadzu LCMS-2020 (Single Quadrupole MS)
 ESI, SIM, positive, Nebulizer gas: 1.5 L/min, Drying gas: 15 L/min
 Heat block temp: 400 °C, DL temp: 250 °C