


**Scherzo SM-C18**  
**Scherzo SS-C18**


150 x 3 mm

Application

**Sodium Halides**  
ハロゲン化ナトリウム類

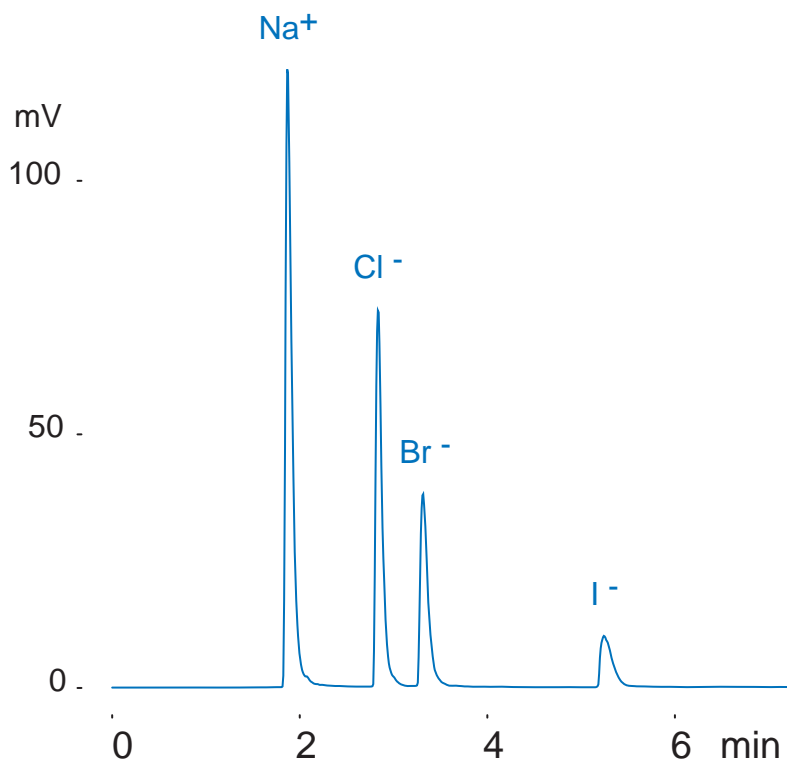

Scherzo SM-C18, 150 x 3 mm  
20 mM ammonium formate  
/ 20 mM formic acid = 95 / 5  
0.4 mL / min ( 12 MPa ), 37 deg.C  
ELSD ( spray chamber 20 deg.C  
drift tube 45 deg.C )  
2.0 uL ( 0.7 ug as sodium salts )

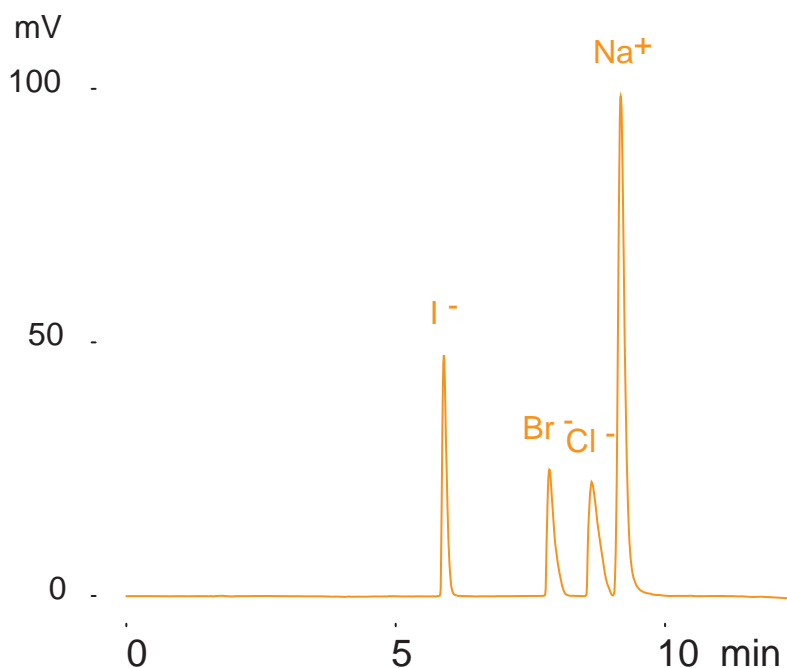

Scherzo SS-C18, 150 x 3 mm  
A: 50 mM ammonium acetate  
/ acetonitrile = 15 / 85  
B: 100 mM ammonium acetate  
/ acetonitrile = 20 / 80  
0 - 100 %B ( 0 - 10 min )  
0.4 mL / min ( 6 MPa ), 37 deg.C  
ELSD ( spray chamber 20 deg.C  
drift tube 45 deg.C )  
4.0 uL ( 1.4 ug as sodium salts )