




Scherzo SS-C18

100 x 2 mm

Application

LC-QTOF/MS: Acyl carnitine

アシルカルニチンのLC-MS/MS分析

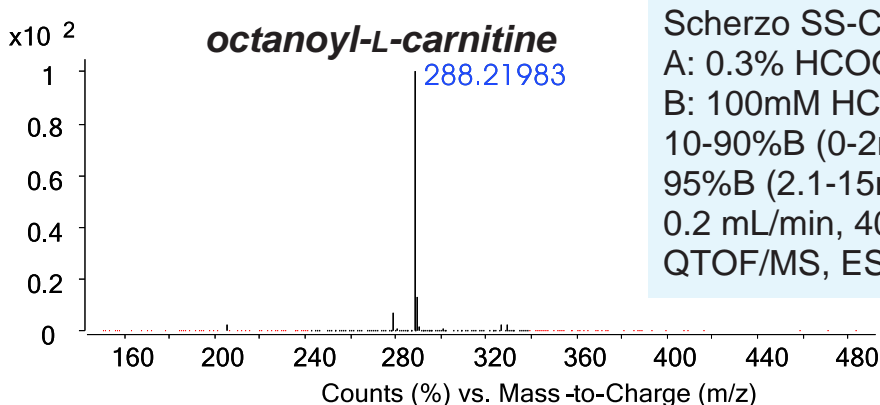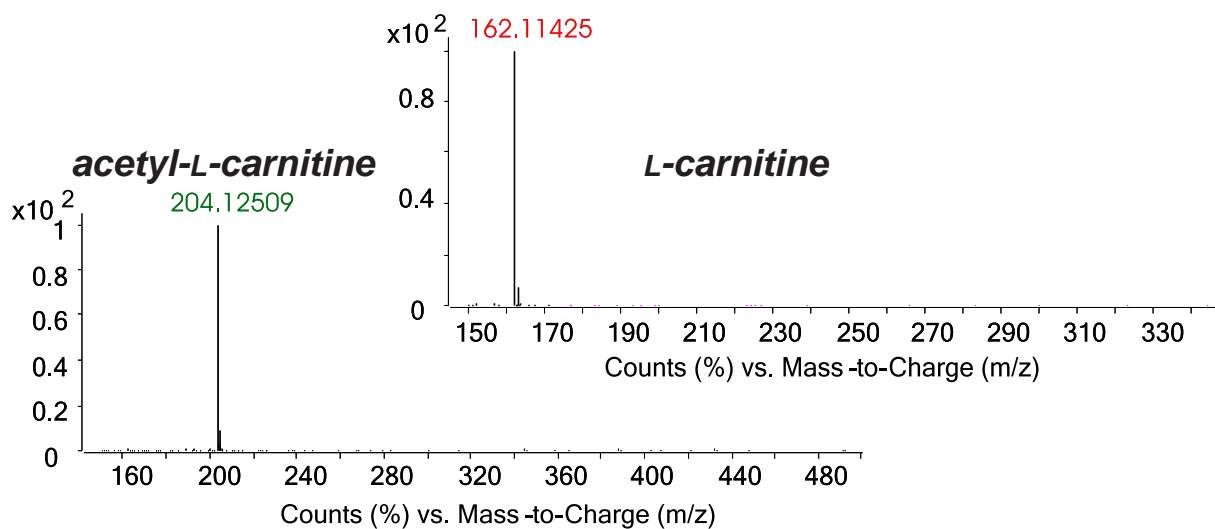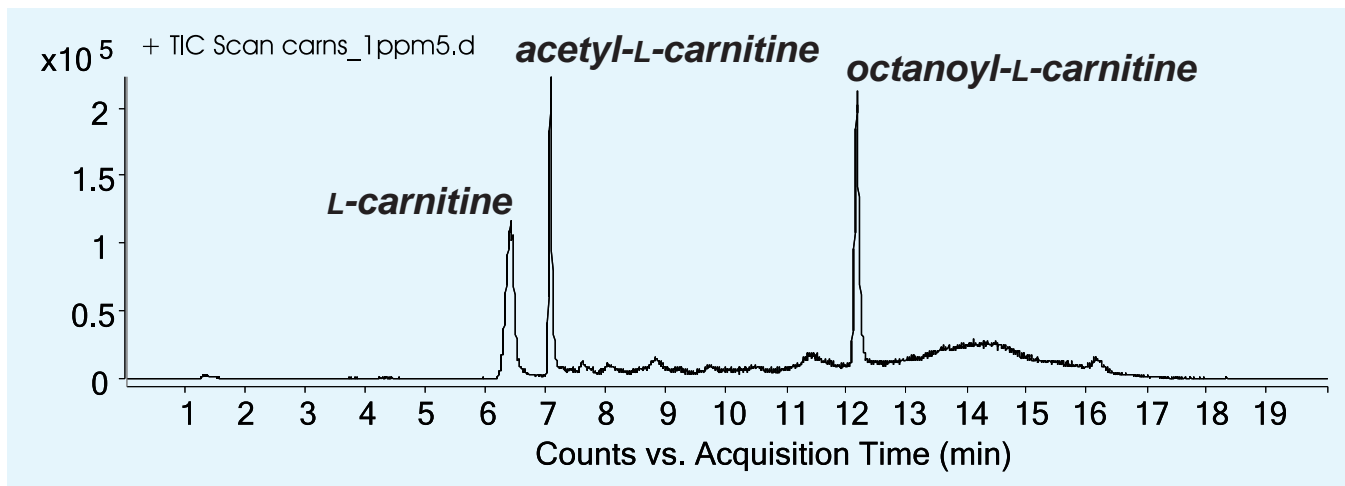

Scherzo SS-C18, 100 x 2 mm
 A: 0.3% HCOOH
 B: 100mM HCOONH₄ / ACN (50 / 50)
 10-90%B (0-2min), 90-95%B (2-2.1min)
 95%B (2.1-15min), 95-10%B (15-16min)
 0.2 mL/min, 40°C, 10 µL injection
 QTOF/MS, ESI, positive

Ref.) I.Nadaoka, E.Hatakeyama, C.Tanada, T.Sakamoto, S.Fukaya, T.Akiba, K.Inoue, Y.Yamano, T.Toyo'oka
 Analytical Sciences, 30, 839 (2014)

Courtesy of Takashi AKIBA, Wakodo Co., Ltd., JAPAN