



Scherzo SS-C18

100 x 3 mm

Application

LC-MS: neurotransmitters (1)

神経伝達物質のLC-MS分析 (1)

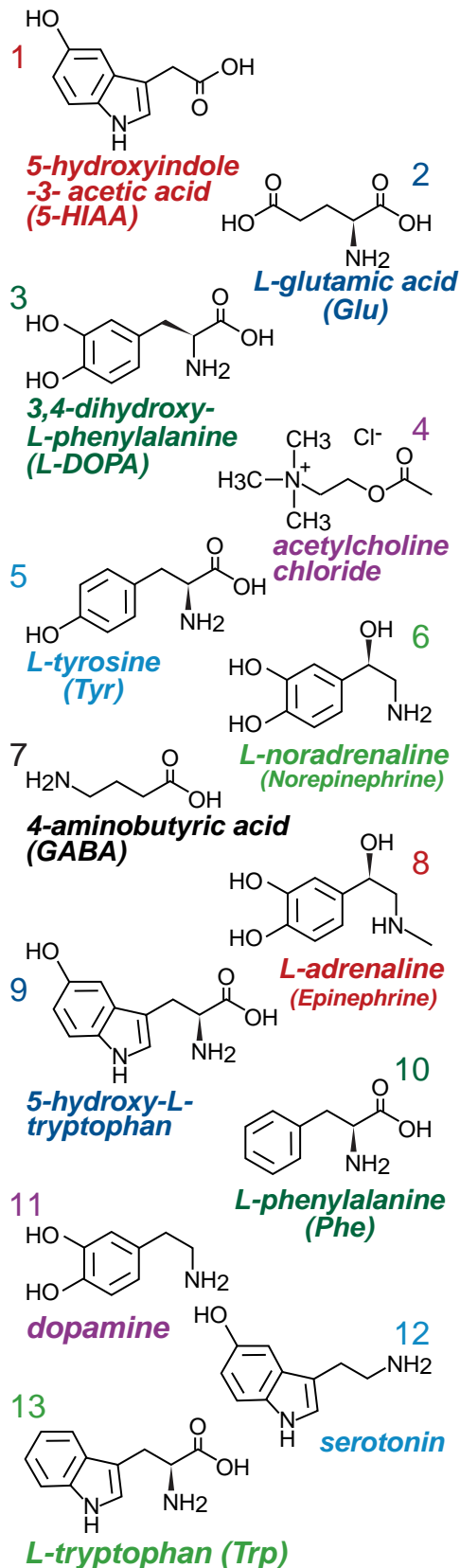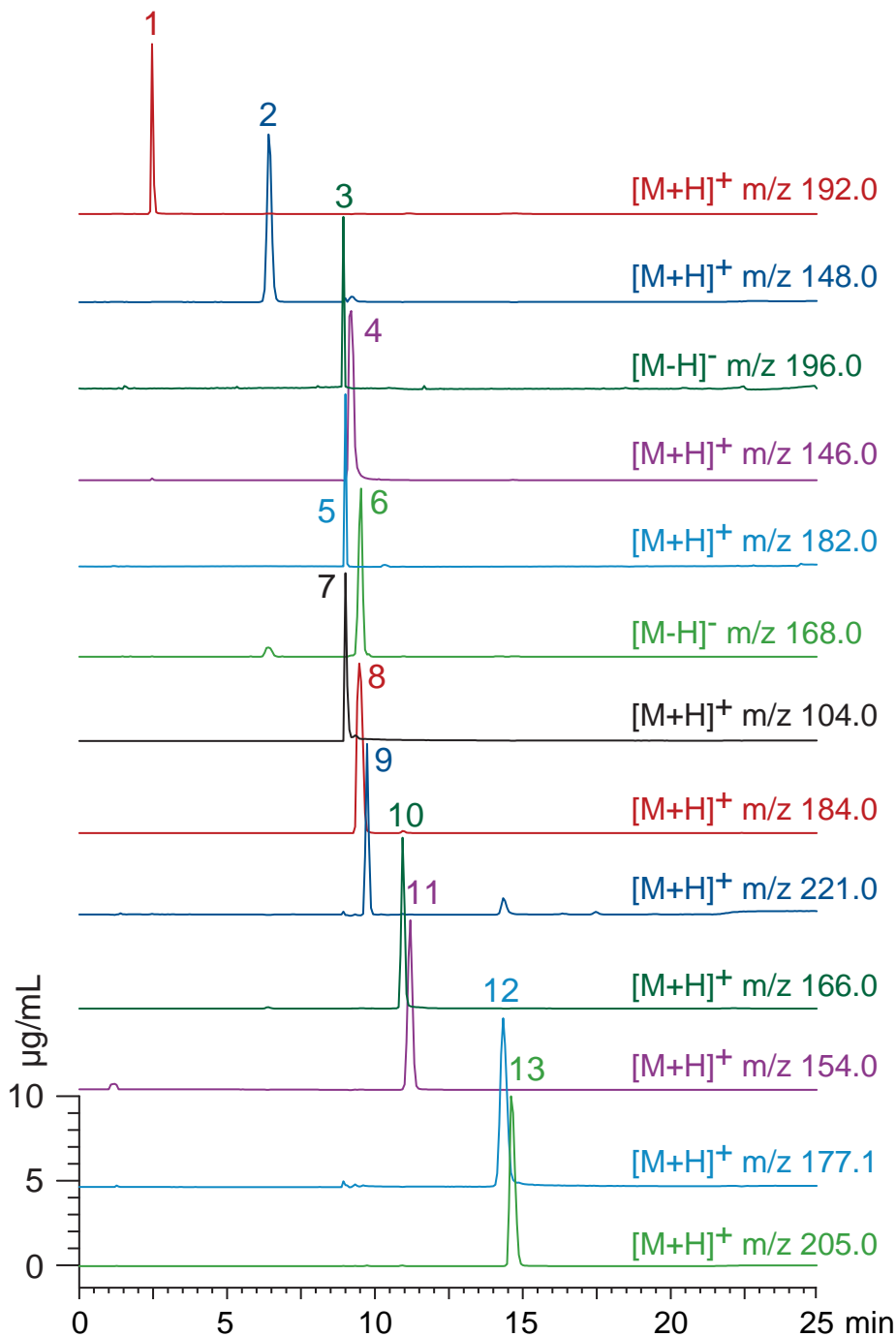

Scherzo SS-C18, 100 x 3 mm

A: methanol / water / formic acid = 35 / 65 / 0.5

B: methanol / 100mM ammonium formate = 40 / 60

0-40 %B (0-20min), 100 %B (20-25min)

0.4 mL/min (9 MPa), 37°C, 20 µL (10µg/mL, 1% HCOOH)

Single Quad. MS (ESI, positive, negative)