






Unison UK-Phenyl
Unison UK-C18

50 x 3 mm

Application

LC-MS: Neuroblastoma related compounds

神経芽細胞腫関連化合物のLC-MS分析

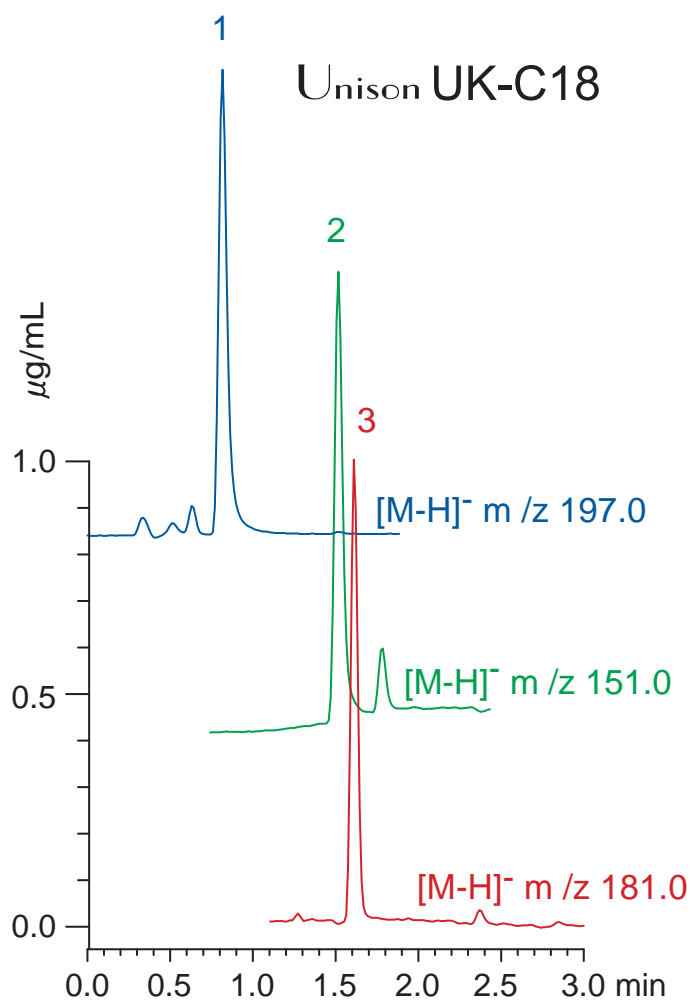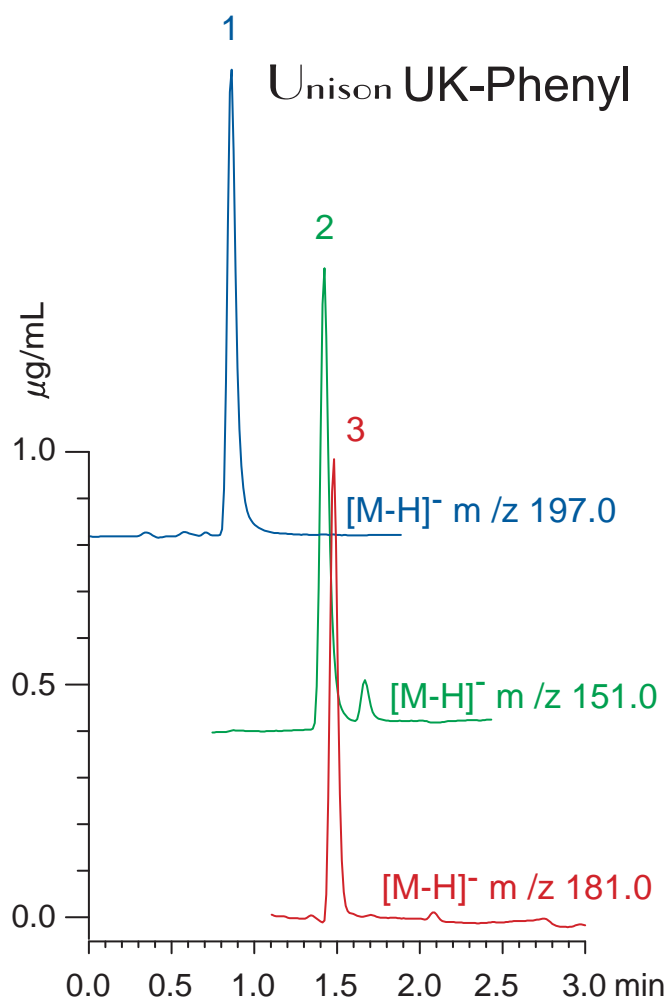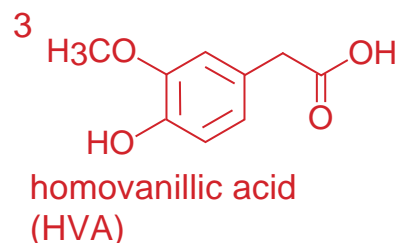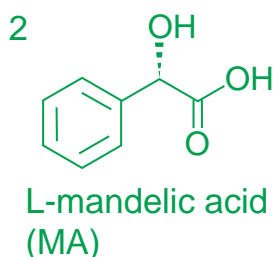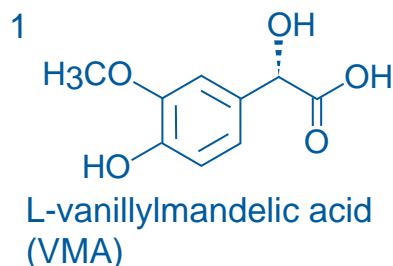

50 x 3 mm

A: water / acetic acid = 100 / 0.5

B: methanol / acetic acid = 100 / 0.5

20-100 %B (0-2 min), 100 %B (2-3 min)

0.5 mL/min (5 MPa), 50°C, 10 μL (1 $\mu\text{g/mL}$, 1% acetic acid)

Single Quad. MS (ESI, negative)