







DACAPO DX-C18

50 x 4.6 mm

Application

Paraben isomers

パラベン異性体

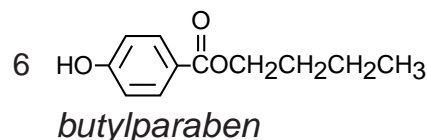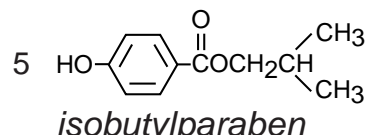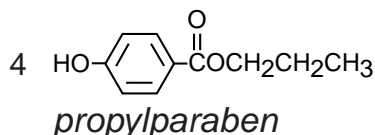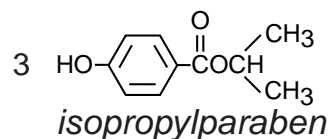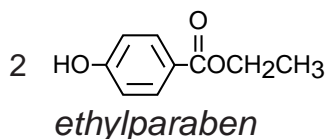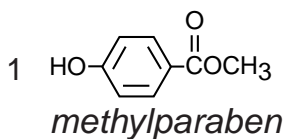

DACAPO DX-C18

Rs = 1.79

5 6

6.6 MPa

50 x 4.6 mm
acetonitrile / water / acetic acid
= 20 / 80 / 0.1
1 mL/min, 37°C, 260 nm
1 μL (0.333mg/mL)

Core-Shell ODS (K)

Rs = 1.32

8.7 MPa

Hybrid Silica ODS (X)

Rs = 1.20

10.3 MPa

0

20

40 min