

**Intrada SEC**


250 x 4.6 mm

Application



**SEC: Glycosaminoglycans (GAGs)**

グリコサミノグリカン(GAGs)のSEC分離



**Glycosaminoglycans (GAGs) Intrada SEC**

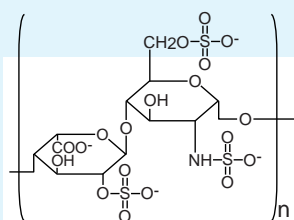


**heparin**

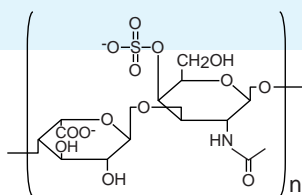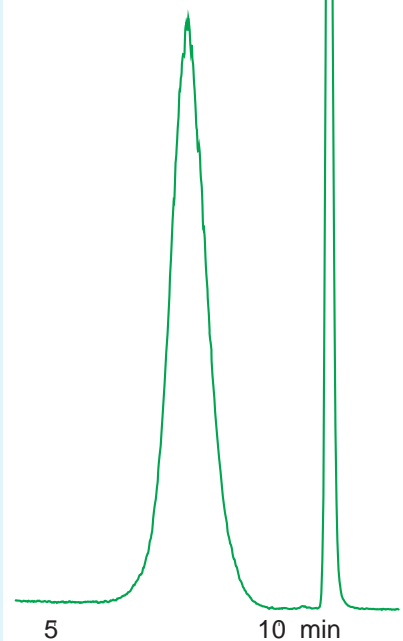

**dermatan sulfate**

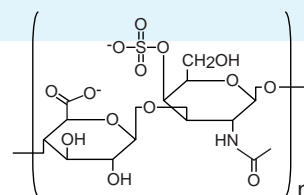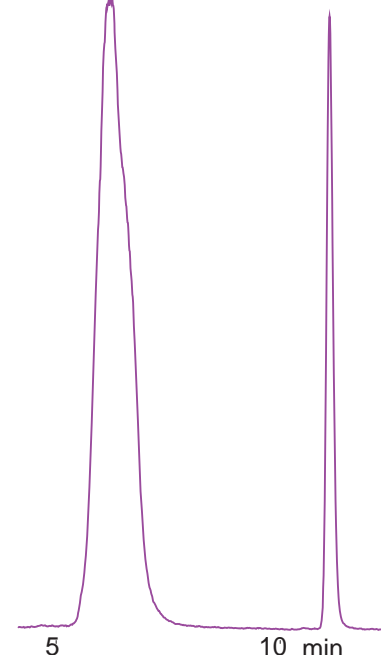

**chondroitin sulfate**

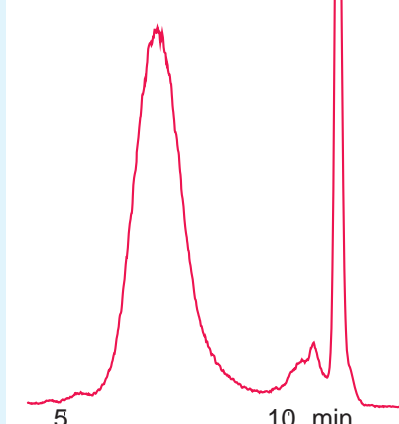

*Intrada SEC, 250x4.6mm  
 100mM ammonium formate  
 0.3 mL/min(5MPa), 37°C  
 ELSD (spray chamber 50°C, drift tube 100°C)  
 10 μL (5mg/mL)*