

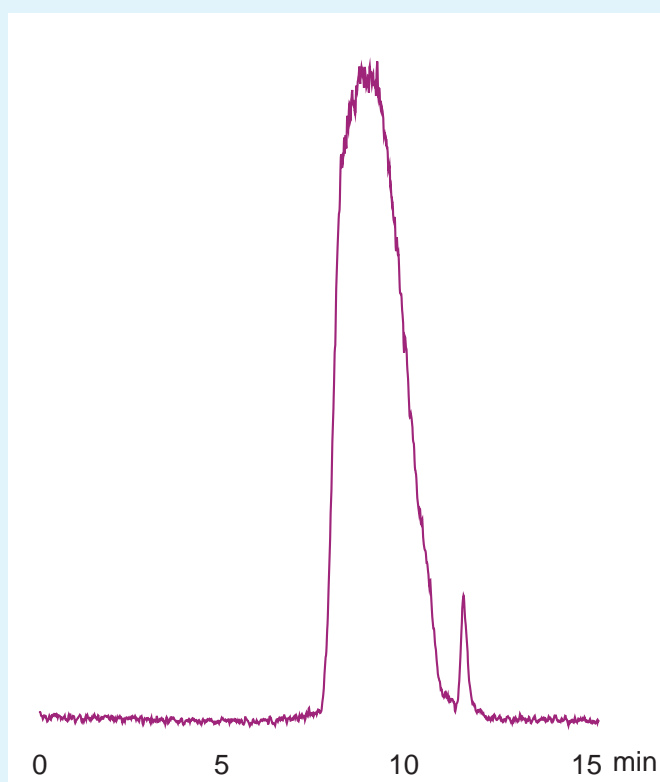
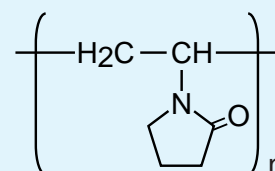
Intrada SEC

250 x 4.6 mm

Application

SEC: Polyvinylpyrrolidone (PVP)

ポリビニルピロリドン(PVP)のSEC分離

**SEC: PVP****Intrada SEC****Polyvinylpyrrolidone K30****MW 40 kDa**

873

Intrada SEC, 250 x 4.6 mm
100mM HCOONH₄ /MeOH = 70 / 30
0.3mL/min (6.3MPa), 50°C
ELSD (spray chamber 50deg.C, drift tube 100deg.C)
10 μL (2mg/mL)