

Intrada SEC

250 x 4.6 mm

Application


SEC: Polyvinyl Acetate (PVAc)

ポリ酢酸ビニル(PVAc)のSEC分離




Polyvinyl Acetate (PVAc)

Intrada SEC

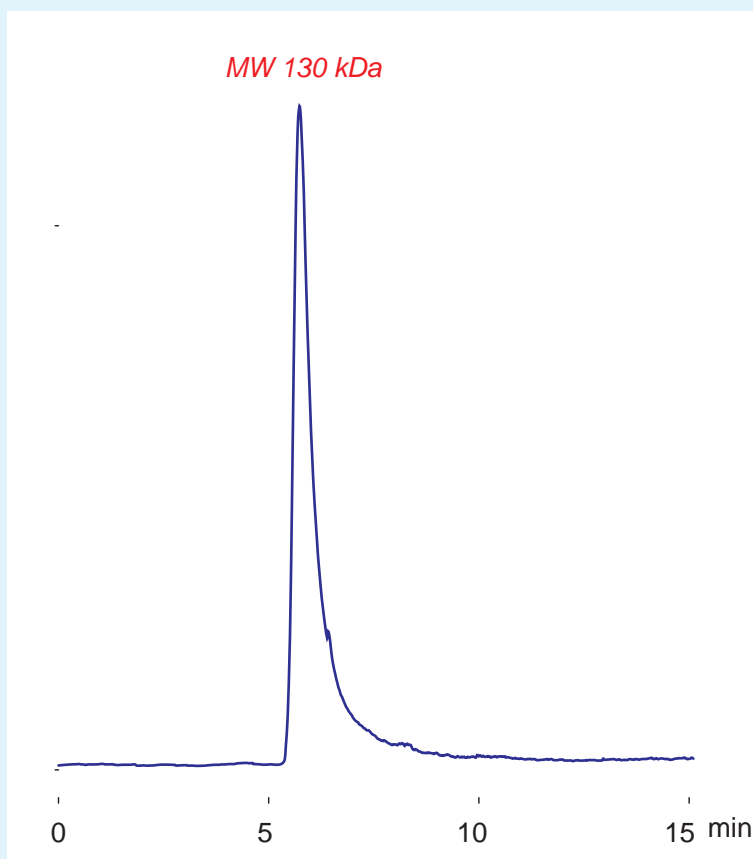

Vinyl Acetate, Polymer

n = 1400 - 1600

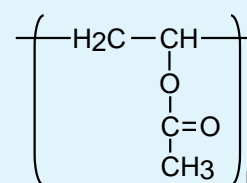

874

Intrada SEC, 250 x 4.6 mm
THF
0.3mL/min (3MPa)
37 °C
ELSD (spray chamber 50°C, drift tube 100°C)
1 µL (2mg/mL)