


Cadenza CD-C18

250 x 4.6 mm


Application

ephedrine, atropine


エフェドリン, アトロピン

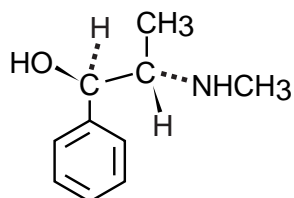

ephedrine

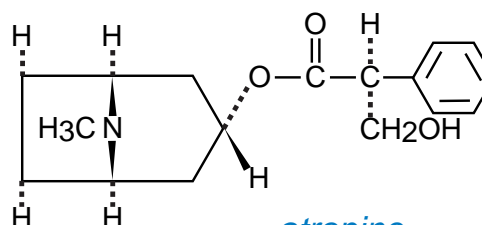

atropine

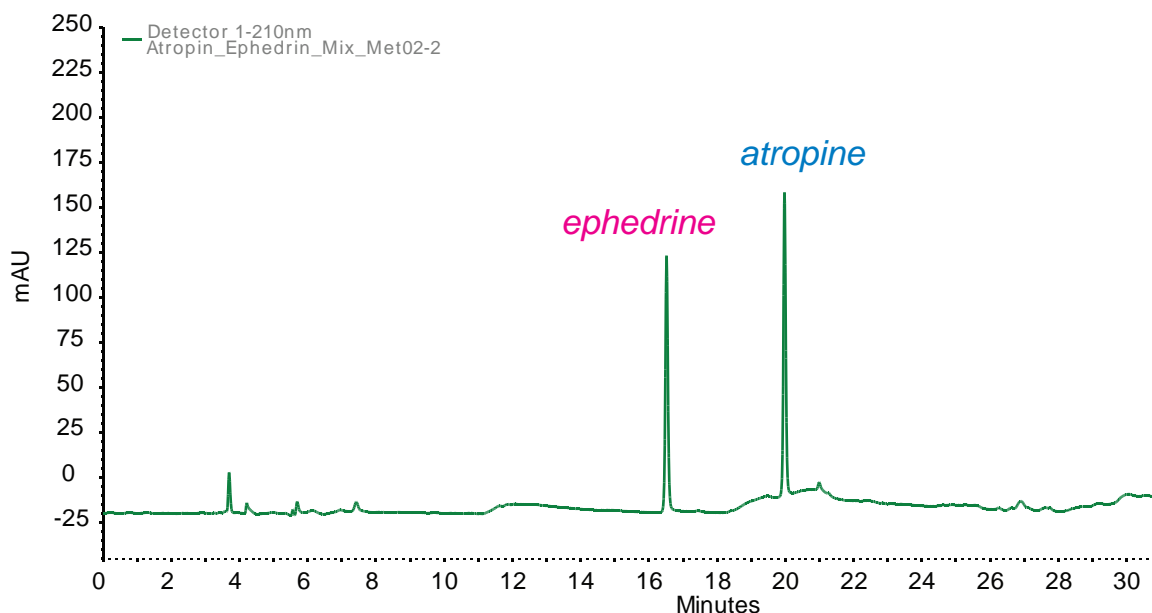

Cadenza CD-C18, 250 x 4.6 mm

A: 0.2% H₃PO₄ in 20mM NaH₂PO₄

B: 0.2% H₃PO₄ in 20mM NaH₂PO₄ / acetonitrile = 20 / 80

5%B (0-5min), 5-70%B (5-30min), 70-5%B (35-35.1min), 5%B(35.1-45min)

0.7 mL/min, 40°C

10 µL (20mg/mL)

210nm

Courtesy of PEAKMAN SP, KOREA