


**Intrada SEC**

250 x 4.6 mm


Application

**SEC: Alkylbenzenesulfonate (LAS)**

陰イオン界面活性剤(アルキルベンゼンスルホン酸)のSEC分離

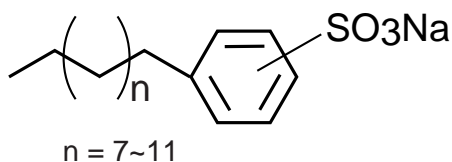

sodium alkylbenzenesulfonates

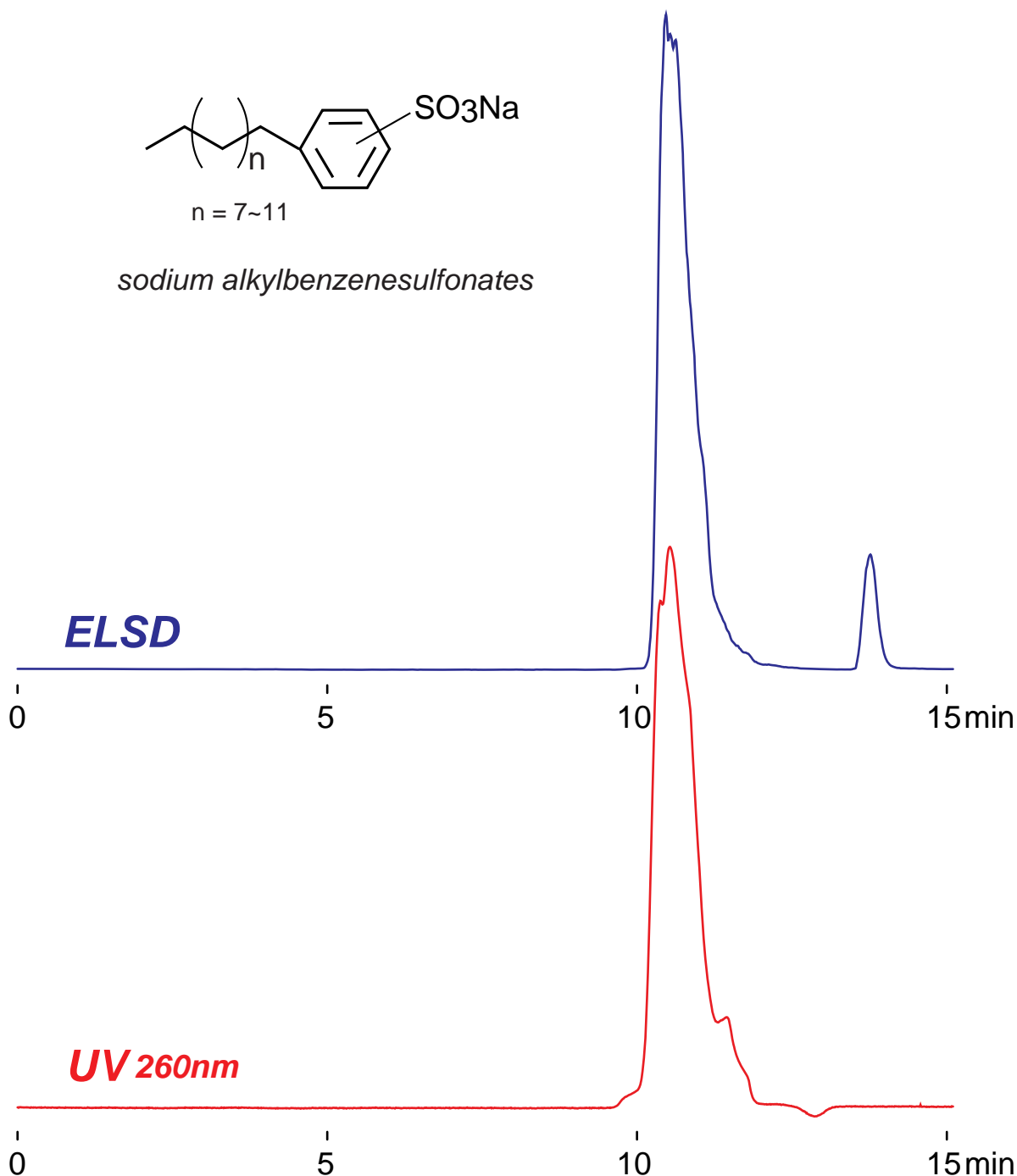

Intrada SEC, 250 x 4.6 mm  
 100mM ammonium formate / acetonitrile = 40 / 60  
 0.3 mL/min (4MPa), 37°C  
 20  $\mu$ L (0.1mg/mL in 100mM HCOONH<sub>4</sub>)