


**Scherzo SM-C18 MF**

100 x 2 mm

Application

**Metal Free Column: picolinic acid**

メタルフリーカラム: ピコリン酸

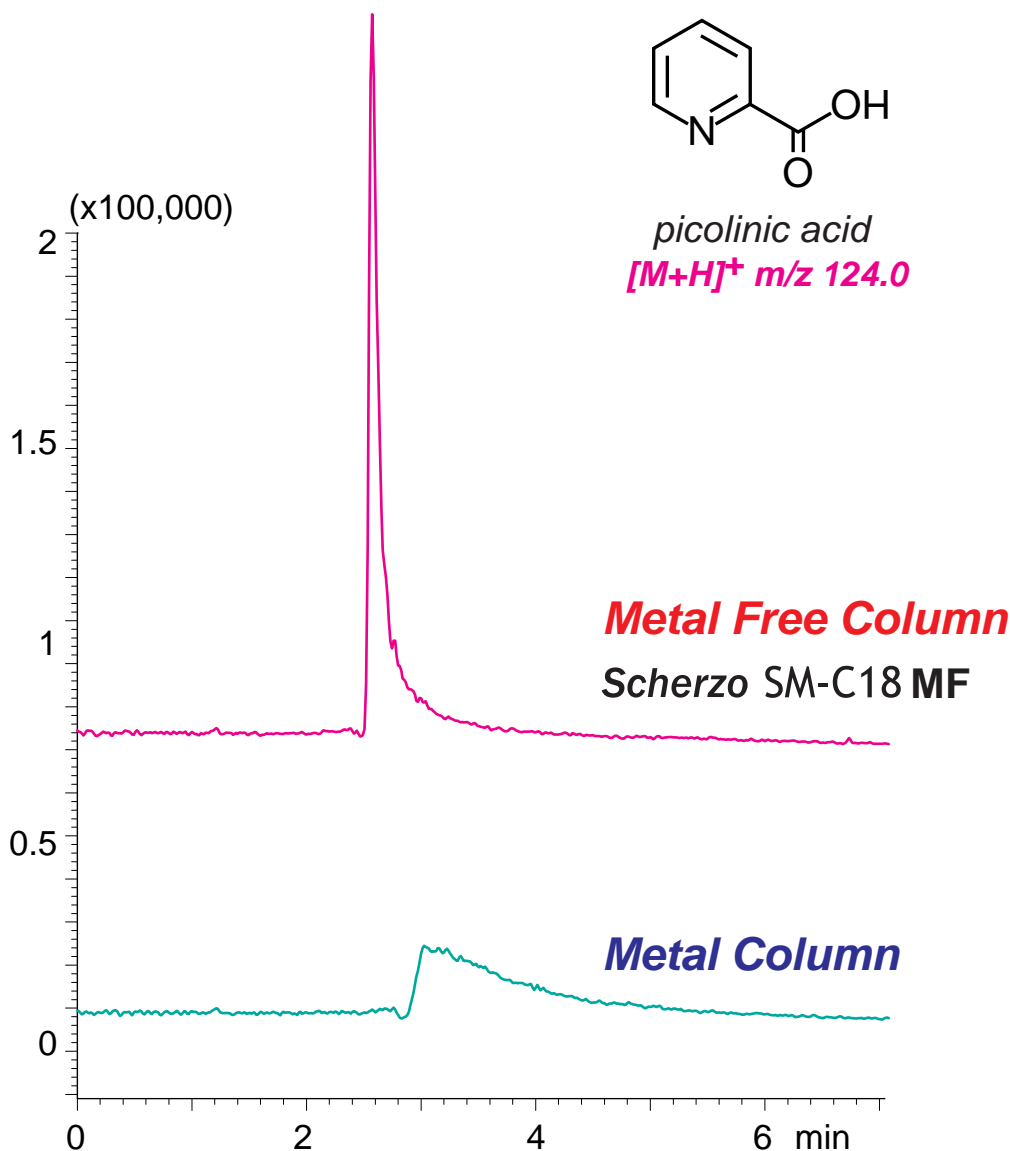

Scherzo SM-C18 MF, 100 x 2 mm  
 A: 50mM ammonium formate  
 B: 100mM ammonium formate / acetonitrile = 50 / 50  
 0-100%B (0-7 min)  
 0.2 mL/min (6MPa), 37°C  
 2 µL (100µg/mL in 1% HCOOH)  
 Single Quad. MS (ESI)