


**Scherzo SM-C18 MF**

100 x 2 mm

Application

**LC-MS: Fe(III)-EDTA and Cu (II)-EDTA**

Fe(III)-EDTAおよびCu(II)-EDTA

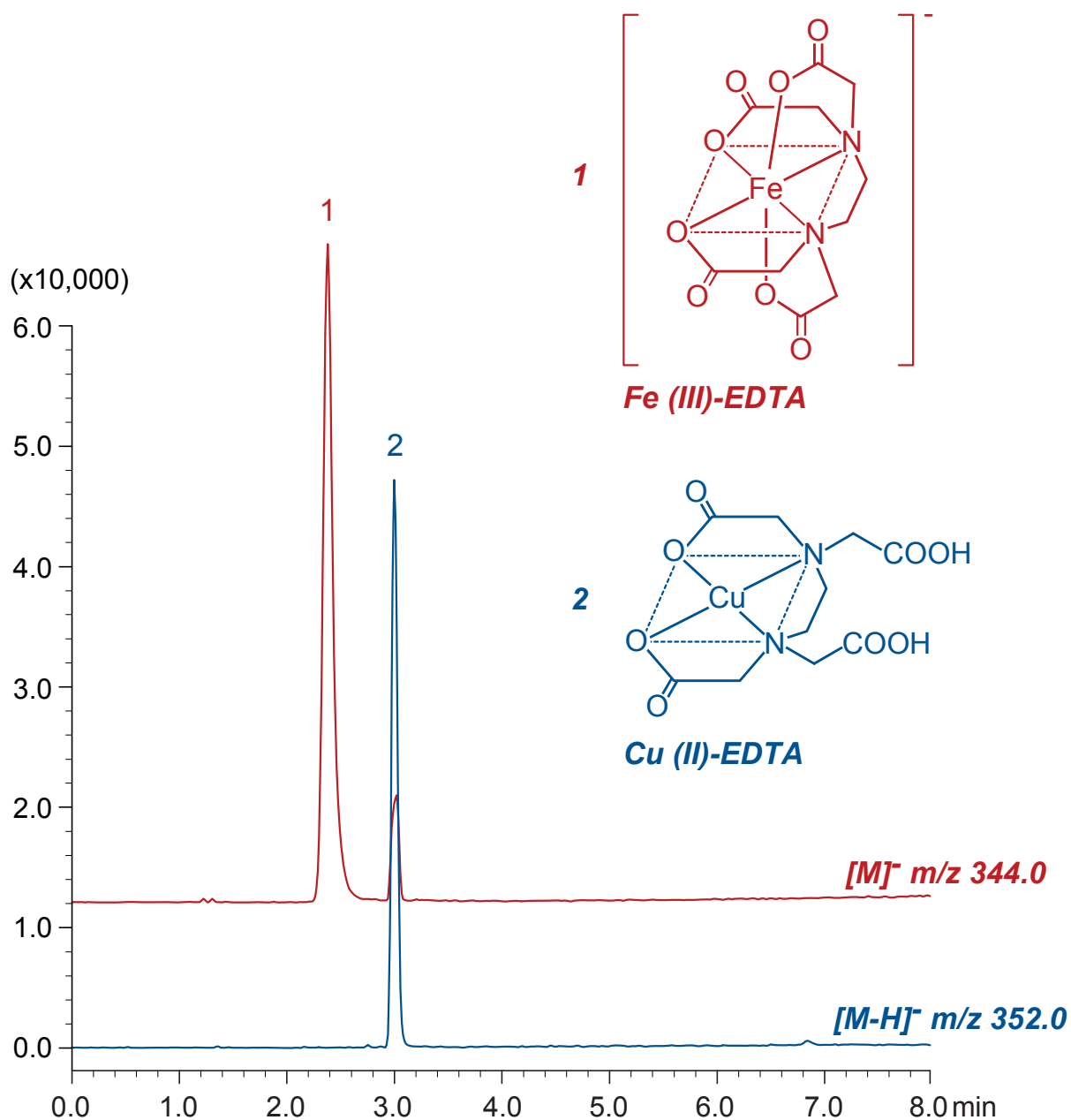

Scherzo SM-C18 MF, 100 x 2 mm

A: 10mM ammonium formate

B: 100mM ammonium formate / acetonitrile = 60 / 40

0-100%B (0-6 min), 100%B (6-8min)

0.2 mL/min (6 MPa), 60°C

2 μL ( 25μg/mL, 5mM ammonium formate), Single Quad. MS (ESI)