

**Cadenza CX-C18**


150 x 3 mm

Application

**Paraben isomers**

パラベン異性体の分離

- 1 methylparaben
- 2 ethylparaben
- 3 isopropylparaben
- 4 propylparaben
- 5 isobutylparaben
- 6 butylparaben

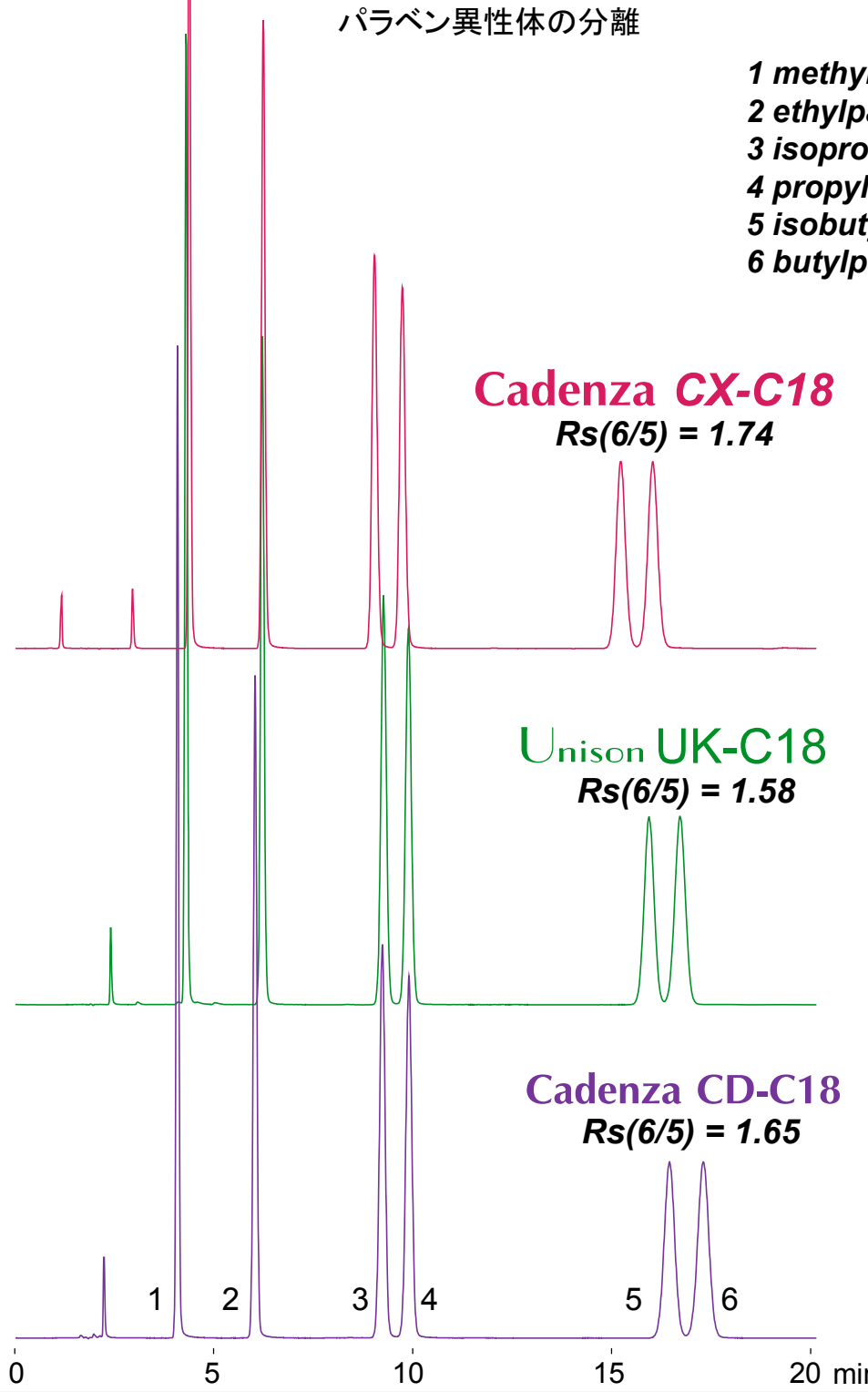

**Cadenza CX-C18**

$R_s(6/5) = 1.74$

**Unison UK-C18**

$R_s(6/5) = 1.58$

**Cadenza CD-C18**

$R_s(6/5) = 1.65$

0 5 10 15 20 min

150 x 3 mm

acetonitrile / water / acetic acid = 40 / 60 / 0.1

0.4 mL/min (8-9MPa), 37°C

1µL (0.33mg/mL), 260nm