






**Scherzo SM-C18 MF**

100 x 2 mm

Application

LC-MS: Organophosphorus pesticides/herbicides

有機リン系農薬のLC-MS分析

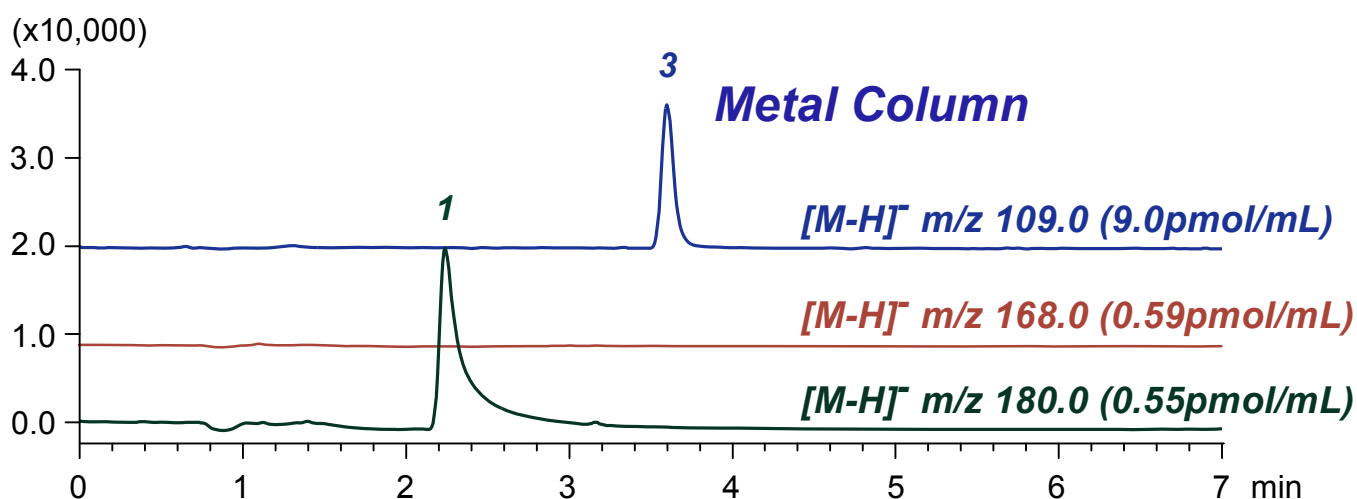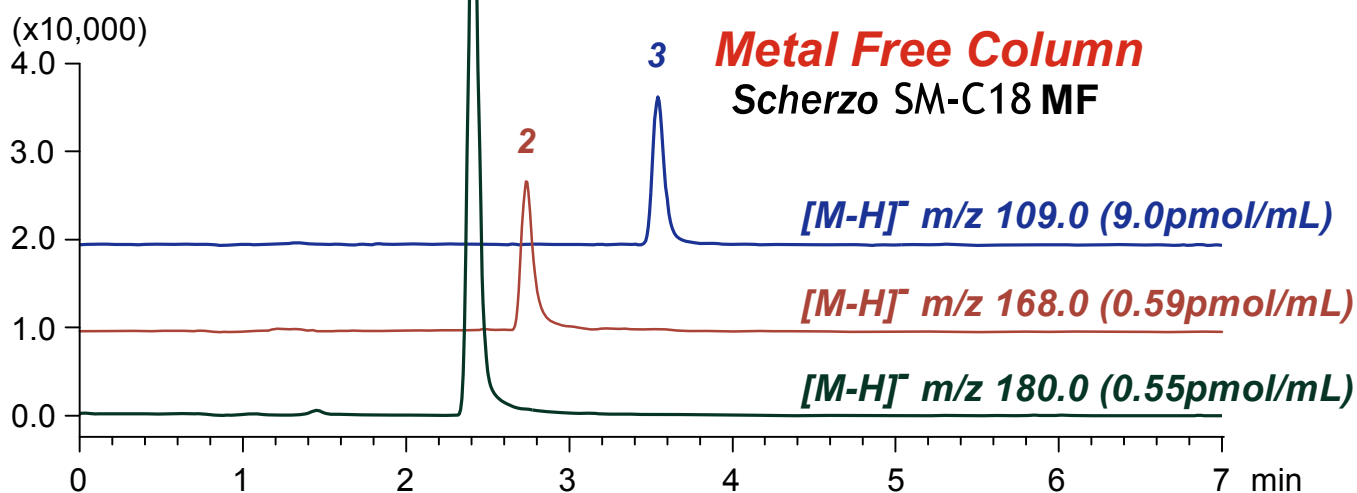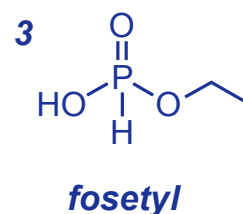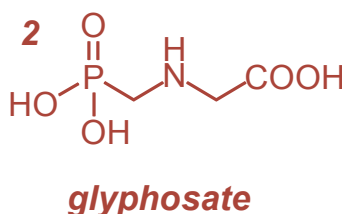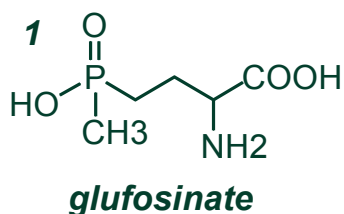

Scherzo SM-C18 MF, 100 x 2 mm

A: methanol / water / formic acid = 20 / 80 / 0.1

B: methanol / 50mM ammonium formate = 20 / 80

0-100%B (0-5min), 100%B (5-7min)

0.2 mL/min (10MPa), 37°C, 20  $\mu$ L (0.55-9.0pmol/mL, water)

Single Quad. MS (ESI)