



Scherzo SM-C18
Cadenza CD-C18

50 x 3 mm

Application

Rizatriptan benzoate (treatment of migraine headaches)

リザトリプタン 安息香酸塩 (片頭痛薬)

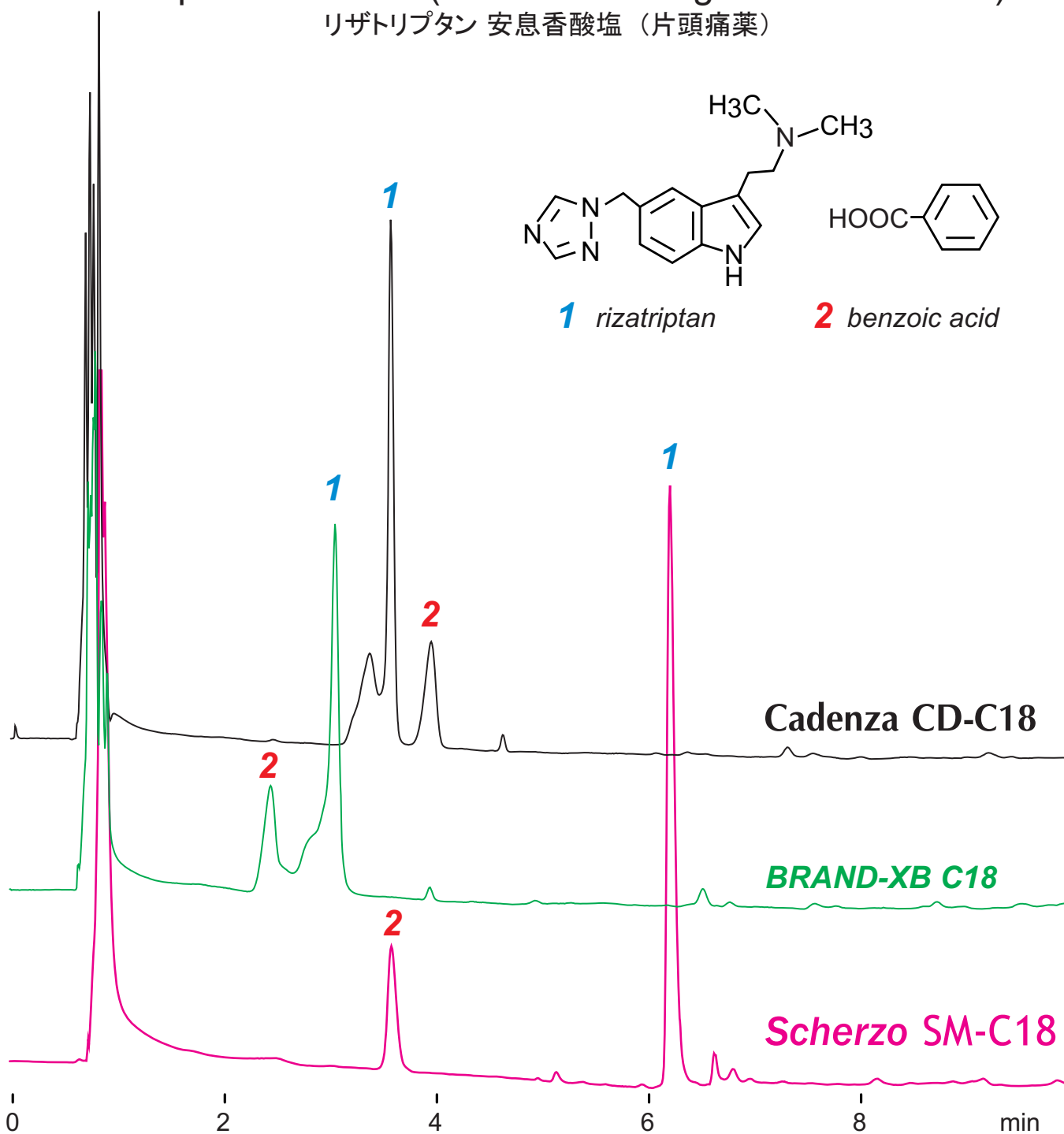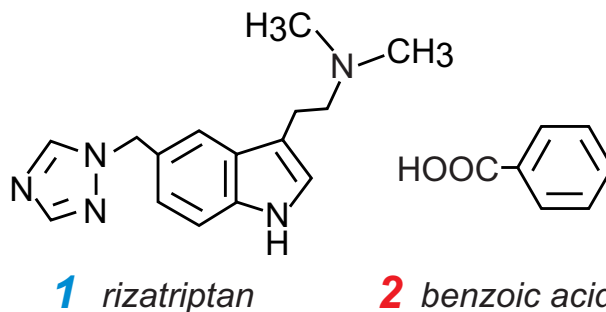

50 x 3 mm
 A: 20mM HCOONH₄
 B: 50mM HCOONH₄ / acetonitrile = 50 / 50
 5%B (0-0.5min), 5-60%B (0.5 -10min)
 0.4mL/min (5MPa), 37°C
 4 μL (0.1mg/mL, DMSO), 260nm