


Cadenza CD-C18
Cadenza CX-C18
Unison UK-C18

75 x 4.6 mm

Application

Tipepidine hibenazate (Antitussive, Expectorant)

チペピジン ヒベンズ酸塩(鎮咳去痰薬)

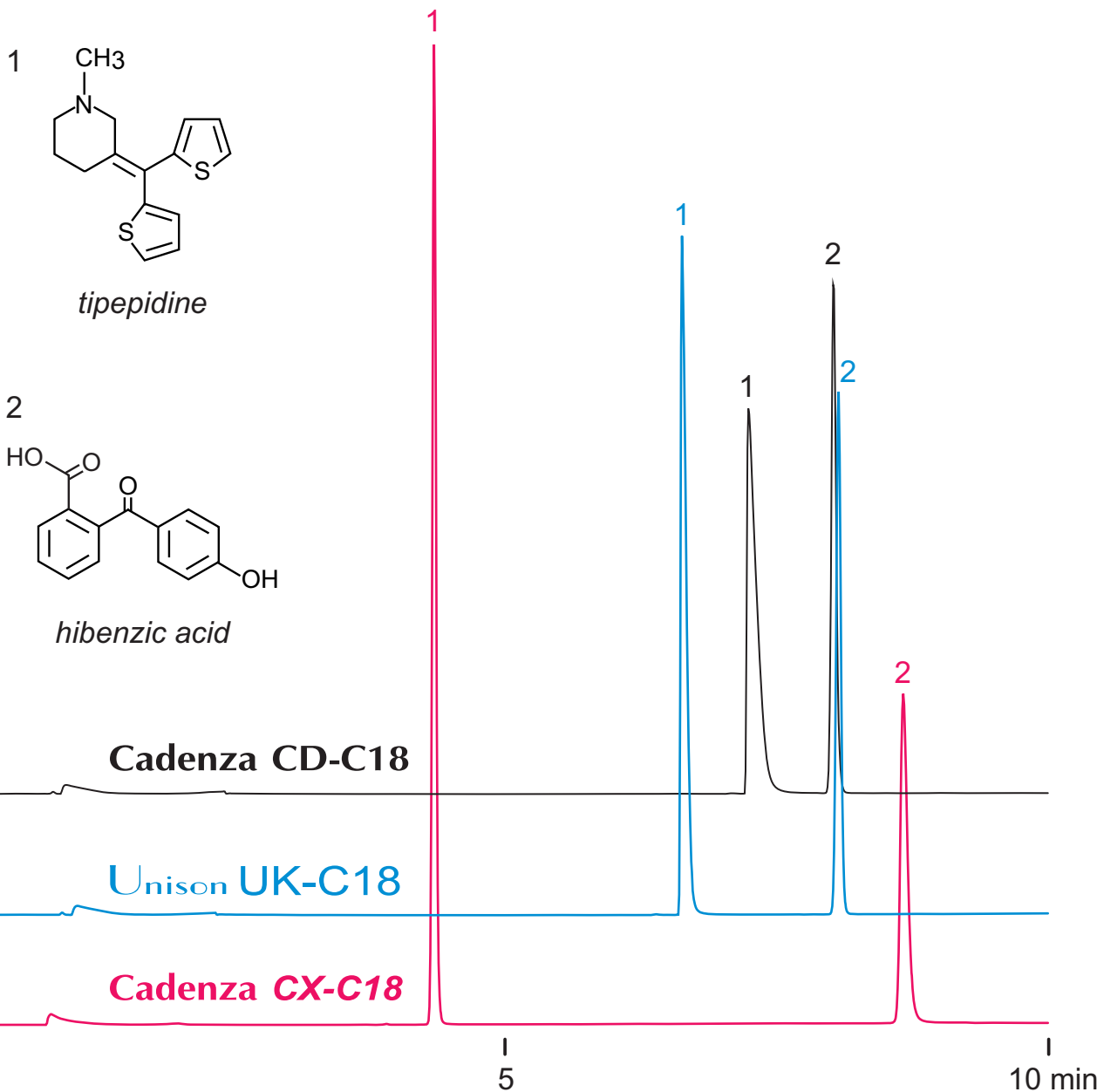

75 x 4.6 mm

A: water / formic acid = 100 / 0.1

B: acetonitrile / formic acid = 100 / 0.1

0-50%B (0-10min)

1.0 mL/min (7MPa), 37°C, 260nm

2 μL (1 mg/mL, DMSO)