


Scherzo SS-C18
 Scherzo SM-C18
 Scherzo SW-C18

50 x 3 mm


Application

Etodolac (nonsteroidal anti-inflammatory drug, NSAIDs)

エトドラク (非ステロイド性鎮痛薬・抗炎症薬)

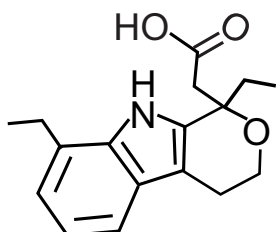

etodolac

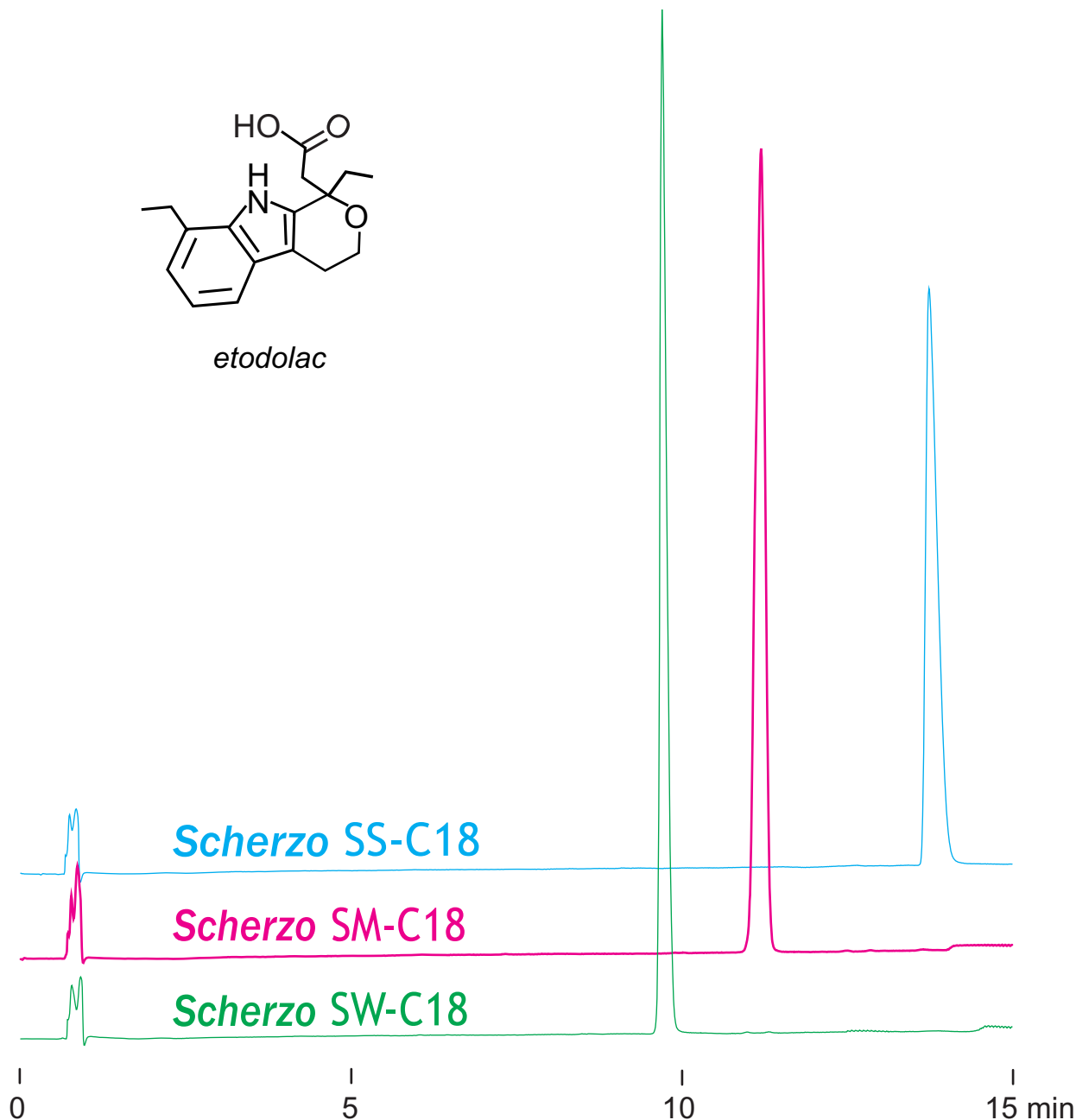

50 x 3 mm

A: 10mM ammonium acetate / acetonitrile = 90 / 10

B: 100mM ammonium acetate / acetonitrile = 50 / 50

0-100%B (0-15 min)

0.4mL/min (5MPa) , 37°C, 260nm

2µL (0.5mg/mL, 0.1% HCOOH in MeOH)