


Scherzo SS-C18
 Scherzo SM-C18
 Scherzo SW-C18

50 x 3 mm

Application

Rizatriptan benzoate (treatment of migraine headaches)

リザトリプタン 安息香酸塩 (片頭痛薬)

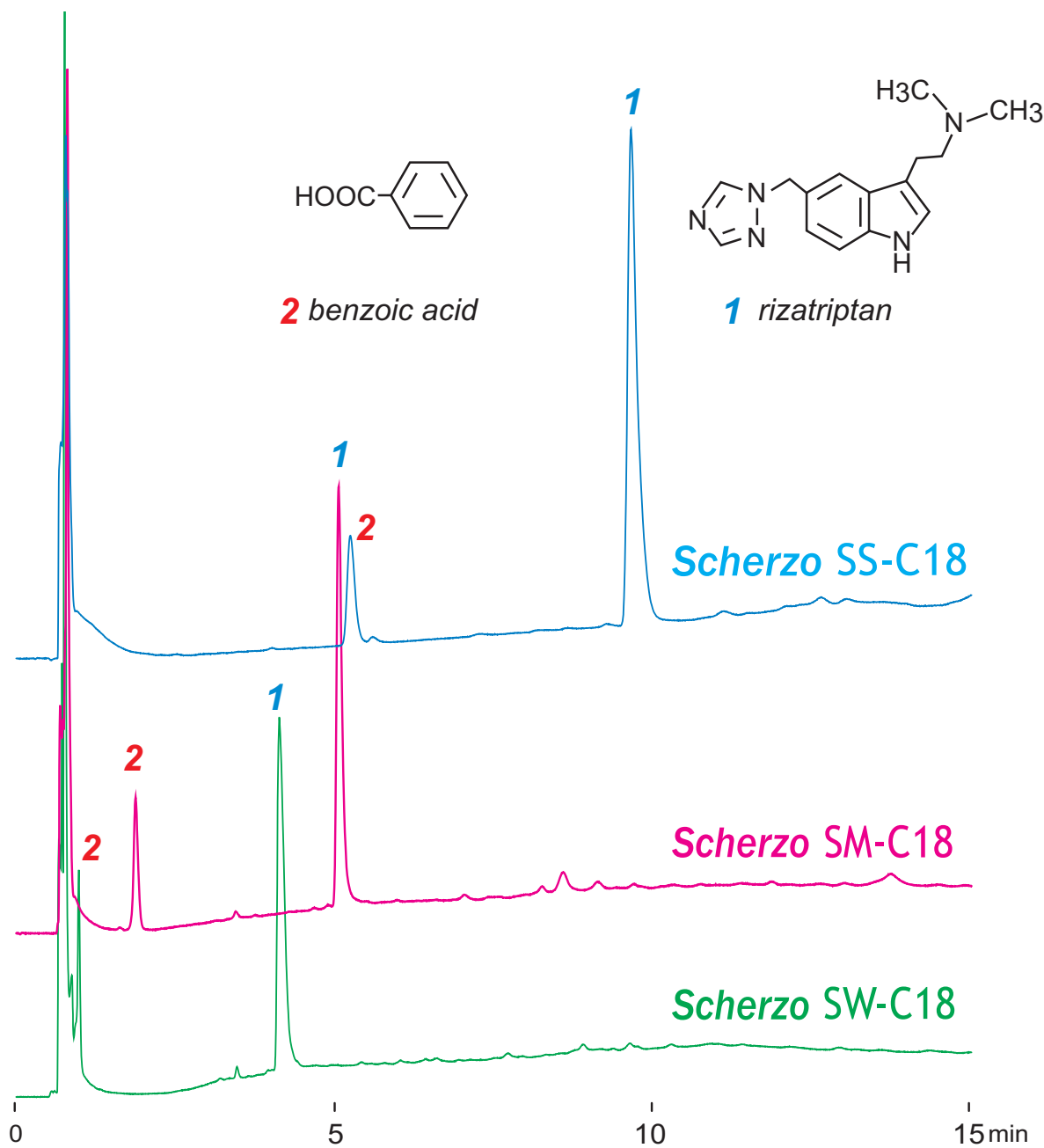

50 x 3 mm

A: 10mM ammonium acetate / acetonitrile = 90 / 10

B: 100mM ammonium acetate / acetonitrile = 50 / 50

0-100%B (0-15 min)

0.4mL/min (5MPa), 37°C, 260nm

2µL (0.1mg/mL, 0.1% HCOOH)