


Scherzo SS-C18  
 Scherzo SM-C18  
 Scherzo SW-C18

100 x 3 mm

Application

LC-MS: Chlorpheniramine maleate (antihistamine)

クロルフェニラミン マレイン酸塩 (抗ヒスタミン薬)のLC-MS分析

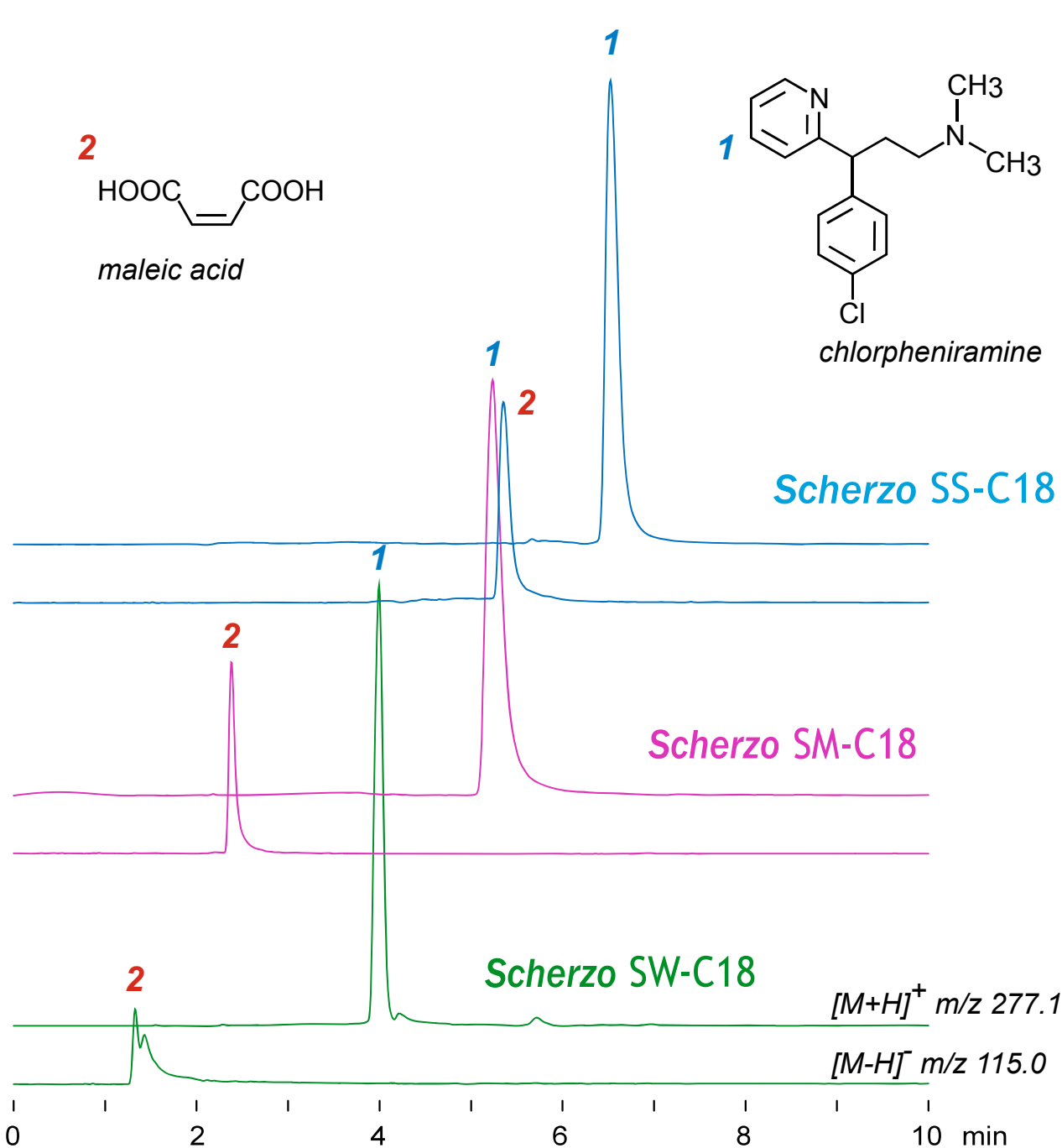

100 x 3 mm

A: 10mM ammonium acetate / acetonitrile = 90 / 10

B: 100mM ammonium acetate / acetonitrile = 30 / 70

0-100%B (0-3 min), 100%B (3-10 min)

0.4mL/min (7MPa), 37°C, MS (ESI, SIM)

5µL (5µg/mL, 50mM AcONH4 / ACN = 70 / 30)