




Scherzo SS-C18  
 Scherzo SM-C18  
 Scherzo SW-C18

50 x 3 mm

Application

**LC-MS: Nelfinavir mesylate (antiretroviral drug)**

ネルフィナビル メシル酸塩 (抗HIV薬) のLC-MS分析

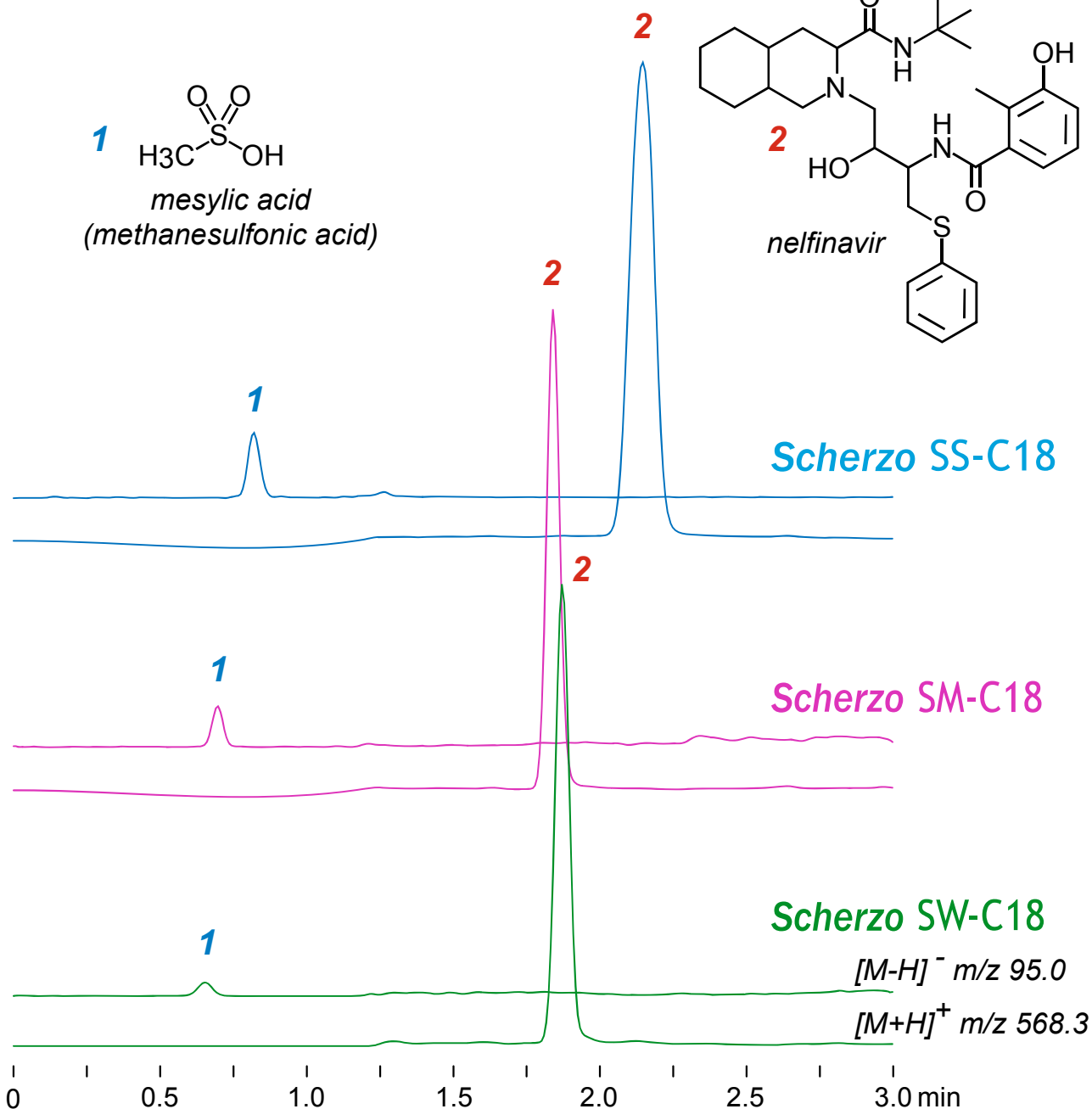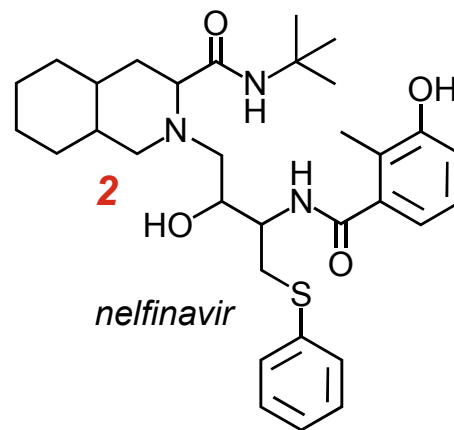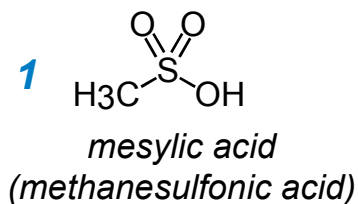

50 x 3 mm

A: 50mM ammonium formate

B: acetonitrile

0-100%B (0-1 min), 100%B (1-3 min)

0.4mL/min (4MPa), 37°C, MS (ESI, SIM)

5µL (2µg/mL, 50mM HCOONH<sub>4</sub> / ACN = 90 / 10)