



Unison UK-C1
 Unison UK-Phenyl
 Unison UK-C8
 Unison UK-C18

50 x 3 mm

Application

LC-MS: Ivermectin (antiparasitic drug)

イベルメクチン(抗寄生虫薬)のLC-MS分析

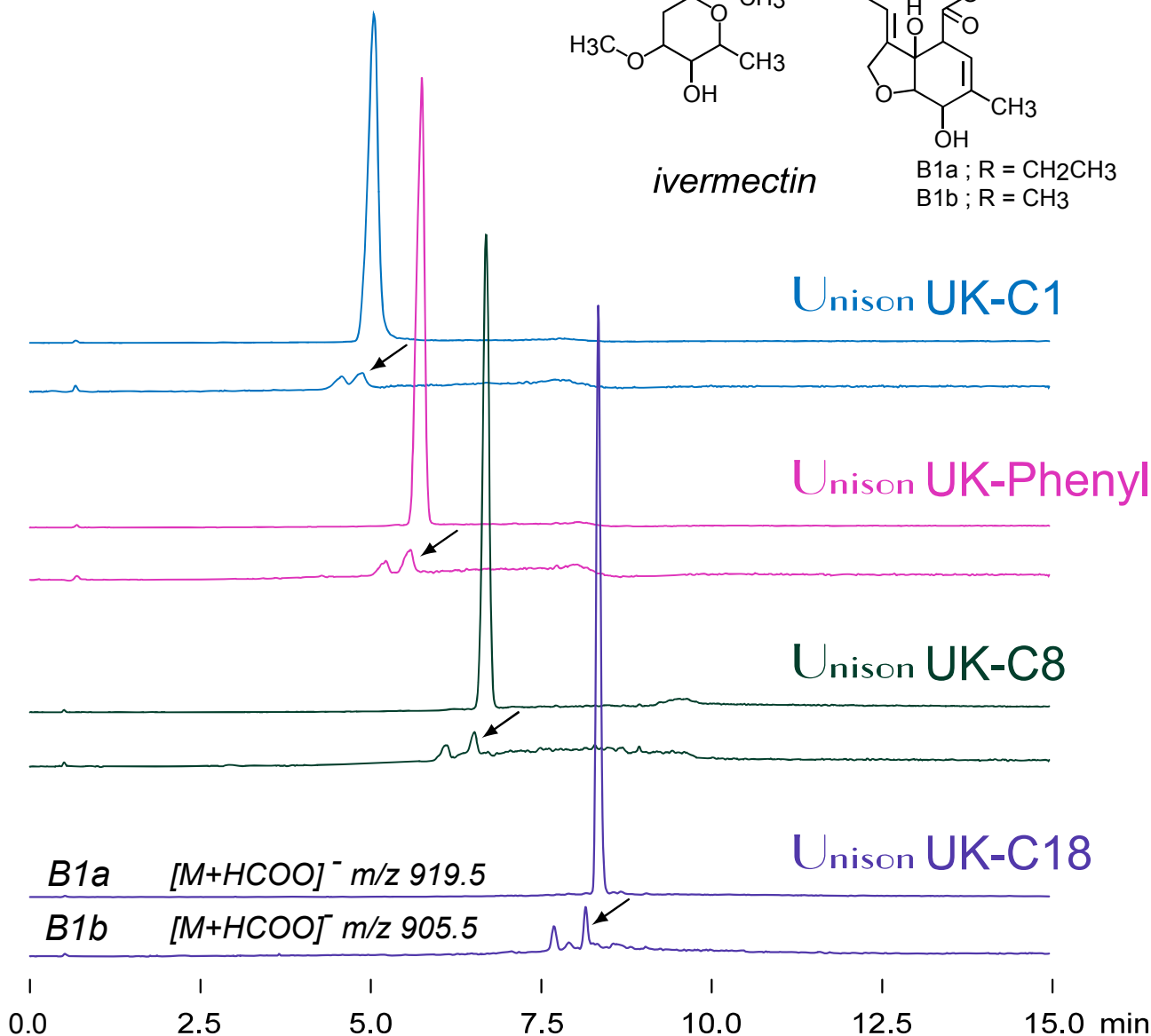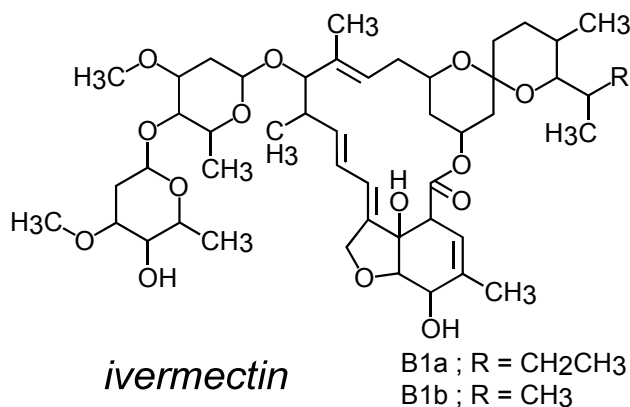

50 x 3 mm

A: water / formic acid = 100 / 0.1

B: methanol / formic acid = 100 / 0.1

50-100%B (0-8 min), 100%B (8-15 min)

0.4mL/min (6MPa) , 37°C, MS (ESI, SIM)

5μL (10μg/mL, EtOH / water / HCOOH = 95 / 5 / 0.1)