


Unison UK-Amino

150 x 3 mm

Application

Favipiravir (RNA polymerase inhibitor)

ファビピラビル (RNAポリメラーゼ阻害剤)

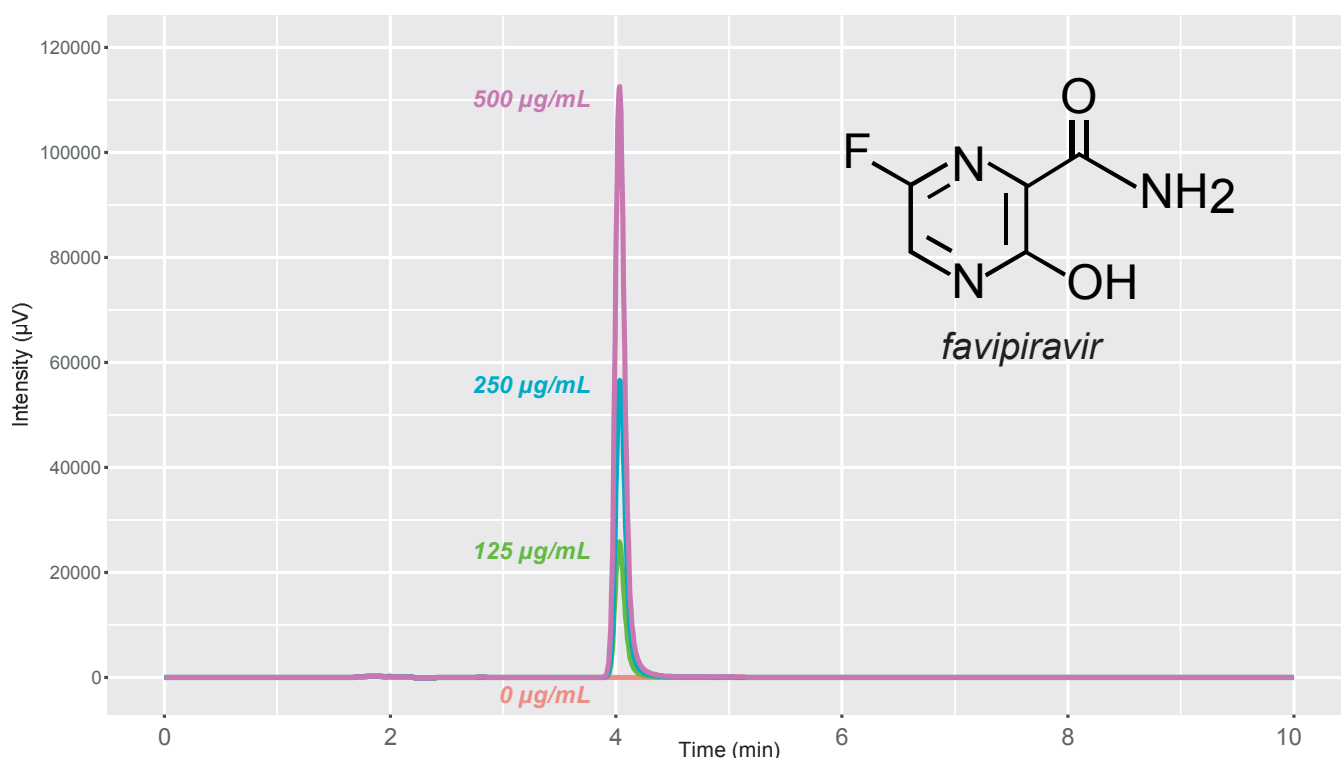

Unison UK-Amino, 150 x 3 mm  
acetonitrile / water / acetic acid = 57 / 43 / 0.3  
0.4mL/min, 60°C, 360nm  
20µL (acetonitrile)

Courtesy of Prof. Dr.TSUJI Yasuhiro and OGAMI Chika, NIHON UNIVERSITY, JAPAN