


Cadenza CX-C18

250 x 4.6 mm


Application

**Rheumatoid arthritis drugs**


関節リウマチ治療薬 (トファシチニブ, フィルゴチニブ)

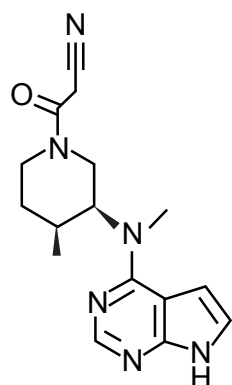

*tofacitinib*

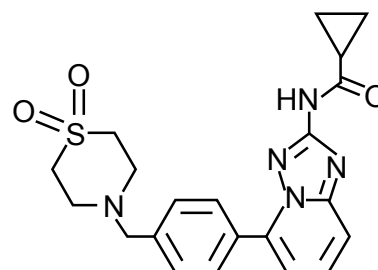

*filgotinib*

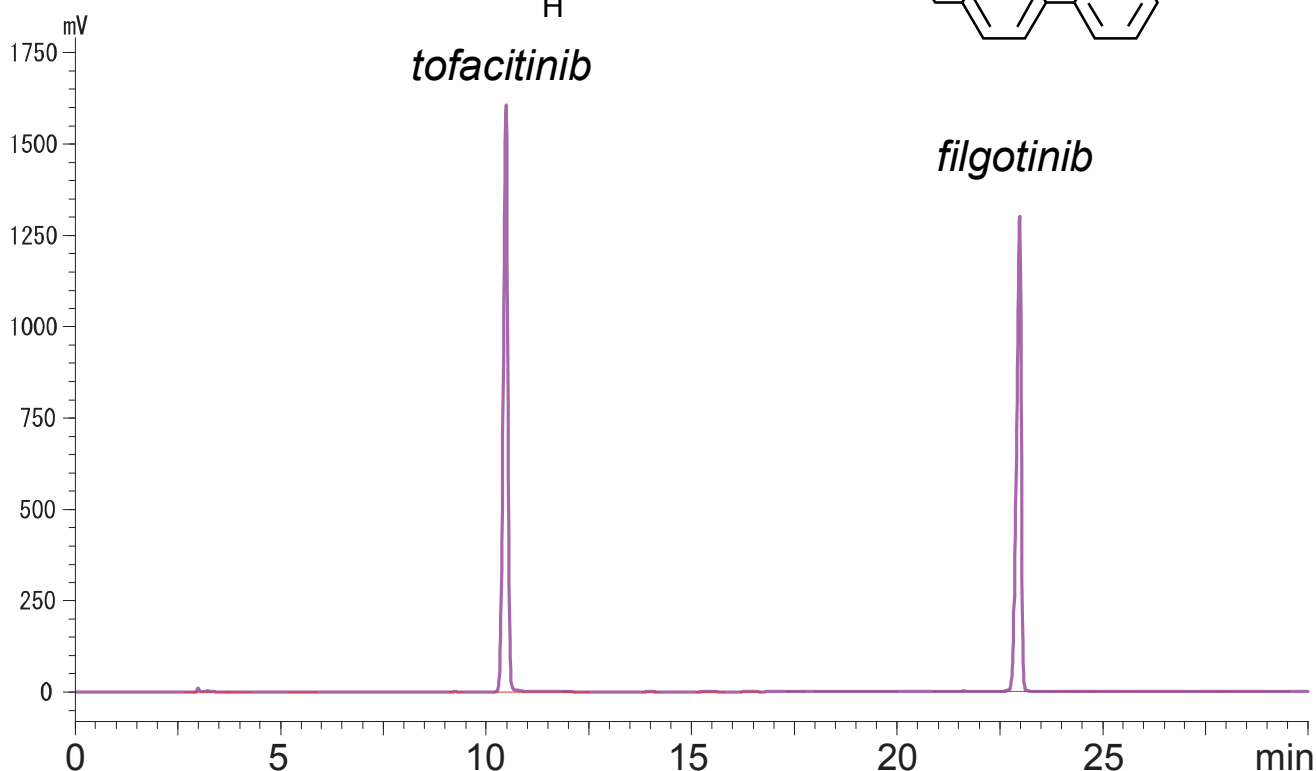

Cadenza CX-C18, 250 x 4.6 mm

A: water / formic acid = 100 / 0.1

B: water / acetonitrile / formic acid = 70 / 30 / 0.1

0-100%B (0-20min), 100%B (20-30min)

1mL/min, 30°C, 287nm

100µL ( 1mg/mL in DMSO, diluted with 1% formic acid )

Courtesy of Prof. Dr. Hideto TO, UNIVERSITY OF TOYAMA (JAPAN)